wau

MORI BLACK

SONIC TOOTHBRUSH

USER GUIDE

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RoHS

Table of contents

Technical Specifications	26
Purpose	26
Device Components.....	27
Operating Procedure	28
Usage Method	28
Cleaning & Maintenance.....	30
Storage.....	30
Charging	31
Disposal.....	31
Warranty.....	33
Questions & Answers	33

Technical Specifications

Parameter	Value
Model	LL-35
Power	3 W
Operating Voltage	5 V
Power Source	Lithium-ion battery
Battery Capacity	750 mAh
Water Resistance Rating	IPX7
Modes & Vibration Frequency	Sensitiv – 31 000 vibrations/min. Gum Care – 35 000 vibrations/min. Polish – 41 000 vibrations/min. White – 45 000 vibrations/min. Clean – 48 000 vibrations/min.
Smart-Timer	2 Minuten
Operating Time	Gentle Mode – 240 Minutes Gum Care – 80 Minutes Polishing – 80 Minutes Whitening – 100 Minutes Cleaning – 200 Minutes
Charging Time	12 hours

Purpose

The MORI BLACK sonic toothbrush by WAU provides effective and gentle cleaning, removing plaque and tartar. The brush operates using high-frequency vibration technology: cleaning is achieved through rapid bristle movement generated by a sonic wave (high-frequency vibration).


Due to the integrated Sensitive Mode, the toothbrush is suitable for sensitive teeth. MORI BLACK offers up to four hours of operation, features five cleaning modes, a memory function, and a zone-change timer for maximum comfort and optimal oral care.

Device Components

1. Brush head
2. Toothbrush body
3. Power on/off/mode switch button
4. Sensitive mode
5. Gum care mode
6. Polish mode
7. White mode
8. Clean mode
9. Low battery indicator
10. Induction charging base

Package contents

- Toothbrush body x1
- Replacement brush heads:
 - with medium-hard bristles x1
 - with soft bristles x1
- Induction charging base x1
- User guidel x1

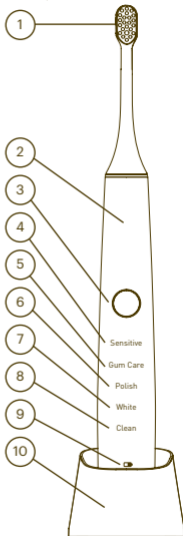Operating Procedure

1. Position the brush head directly above the motor drive shaft and attach it.

Important: Do not screw the brush head on during installation.


2. To turn on the toothbrush, briefly press the power/mode button.

3. Each subsequent short press switches the modes in the following order: Sensitive, Gum Care, Polish, White, Clean.

4. Every 30 seconds, the toothbrush pauses briefly to remind you to change the brushing area.

5. After 2 minutes, the automatic operating cycle will end and the toothbrush will switch off.

6. To turn off the toothbrush manually, press the power/mode button.

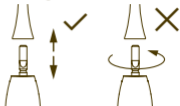Usage Method

1. Attach the brush head.

2. Moisten the brush head with water and apply the required amount of toothpaste.

3. Turn on the toothbrush.

4. Place the brush head in your mouth. Hold it at a 45-degree angle to the gums. Press the power button and select a comfortable operating mode.


5. Begin brushing. Move the brush head along the teeth using sweeping motions from the gum line toward the edge of the tooth. In vier Quadranten unterteilt werden — putzen Sie jeden Quadranten 30 Sekunden lang.

6. Dentists recommend brushing for at least 2 minutes. The oral cavity can be divided into four quadrants — brush each quadrant for 30 seconds.

To support this recommendation, the toothbrush includes a 30-second interval timer: it briefly pauses to signal the need to switch to the next area.

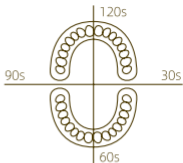After 2 minutes, the automatic operating cycle ends and the toothbrush switches off.

7. After brushing, rinse the brush head and allow it to air dry. The toothbrush body can also be rinsed with water.

Features of the MORI BLACK sonic toothbrush

Brush Head Replacement Indicator

It is recommended to replace the brush head every three months. The medium-soft bristles of the original MORI electric toothbrush head gradually fade in a heart-shaped pattern after three months of use, reminding you that it is time to replace the head.

Memory Function

The toothbrush remembers the last mode you used for more than 5 seconds. The next time you turn it on, it will start in that same mode.

Single-Button Control

Simple and convenient: one button is used for switching modes and turning the toothbrush on and off.

Cleaning & Maintenance

— Brush Heads for the MORI BLACK Sonic Toothbrush

Rinse the brush head with water after each use and allow it to air dry.

To extend the service life of the product, use replacement brush heads specifically designed for the MORI BLACK sonic toothbrush (model LL-35).

— Toothbrush Body

The toothbrush body may be cleaned with water; afterward, allow the device to air dry.

Do not use water or acidic cleaners at temperatures above 40°C.

Storage

— Do not drop the device

Discontinue use if the brush body or brush head is damaged.

— Keep away from heat sources

Do not use the charging base in areas with elevated temperatures.

— Do not assemble or disassemble the MORI BLACK sonic toothbrush yourself

If the toothbrush requires repair, please contact a service center authorized by the manufacturer to handle consumer claims.

Charging

- The MORI BLACK sonic toothbrush is pre-charged before leaving the factory. Please recharge the toothbrush once the battery is fully depleted.
- The low battery indicator flashes red when the device needs charging. (The icon on the toothbrush body does not display the battery level; it only signals a critically low charge).
- While charging, the MORI BLACK sonic toothbrush body will pulse gently (slowly lighting up and dimming). When fully charged, the pulsing stops.
- Do not place metal objects on the charging base during charging to avoid accidents.
- Please recharge the toothbrush only after the battery is fully discharged. This helps extend the lifespan of the lithium-ion battery.

Disposal

Do not dispose of this product with household waste. Please protect the environment and take the device to an appropriate recycling center.

Disposal of the Lithium-Ion Battery

The built-in rechargeable battery contains substances that can harm the environment.

When disposing of the MORI BLACK sonic toothbrush, always remove the battery and dispose of it at a designated battery recycling facility.

If you experience difficulties removing the battery, you may contact a service center authorized by the manufacturer to assist with this procedure.

Remove only the battery when disposing of the toothbrush. Ensure the battery is fully discharged before removal.

Procedure:

1. Disconnect the toothbrush from any power source.
2. Stop using the device.
3. Remove the bottom cover of the toothbrush.
4. Take out the internal mechanism.
5. Remove the battery.

Important: After removing the battery, do not reconnect the toothbrush to a power source.

Names and Content of Hazardous Substances in the Product

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This table is prepared in accordance with the SJ/T 11364 standard.

○ indicates that the content of the hazardous substance in all homogeneous materials of the component is below the limit requirements specified in the GB/T 26572 standard.

×

× indicates that the content of the hazardous substance in at least one homogeneous material of the component exceeds the limit requirements of the GB/T 26572 standard.

This mark identifies the safety of using the product for the environment for a period of 10 years. In this risk warning, the product contains toxic and hazardous substances within normal limits that may cause environmental pollution or mutation. Using this product within the specified period will not pollute the environment or harm the user or others. When disposing of the product after the expiration of this safety period, local regulations and laws for waste recycling must be followed.

Garantie

The MORI BLACK sonic toothbrush by WAU is covered by a 1-year warranty. During the warranty period, the device is eligible for replacement or repair in the event of any technical defects. This warranty does not apply to malfunctions caused by improper use, transportation damage, or unauthorized repairs performed by the user.

Questions and Answers

Q: The brush vibration power becomes weaker and weaker.

A: Charge the device fully or replace the charger. Charge on time.

Q: After charging, the indicator continues to show low charge.

A: Make sure the electric toothbrush was kept vertical during charging, which may affect the quality of the charge. Also, it is possible that the battery life has expired. Contact the authorized manufacturer for consumer claims.

Q: The electric toothbrush was not used for a long time before the first use. When trying to turn it on, the body does not light up, as it should during charging.

A: If the battery has not been used for more than 6 months, start the charging process. Use the original charger, suitable for the standard charger input voltage.

Q: What if the brush still does not turn on?

A: Check if it is fully charged; it may still have a low battery level. Try to turn it on after charging again. If the problem persists, contact the authorized manufacturer for consumer claims.

MORI BLACK

SONIC TOOTHBRUSH

Model: **LL-35**

Material: **PC plastic**

Dimensions (H×W×D): **258×55×47 mm**

Weight: **200 g**

Voltage: **4 V**

Power: **3 W**

Origin: **PRC**

Authorized manufacturer representative: PERFECT SQUARE LIMITED 54/2 Fairclough Rd, Beach Haven, Auckland, Auckland, 0626, New Zealand

e-mail: info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com

RoHS