wau

NEO

PLASMA DEVICE

USER GUIDE

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Table of Contents

Technical Specifications	13
Purpose.....	14
Device Components.....	14
Operating Procedure.....	15
Usage Method	16
Cleaning & Maintenance.....	17
Storage.....	17
Charging.....	17
Precautions	18
Warranty.....	18


Technical Specifications

Parameter	Value
Model	BPW-002
Materials	ABS plastic, PC plastic, silicone
Dimensions (HxWxD)	Device with base: 189x50x45 mm
Weight	Device with base: 183 g
Power source	Lithium-ion battery
Battery capacity	1000 mAh
Operating voltage	3.7 V
Power	5 W
Charging time	3 hours
Operating time (from full charge)	60 minutes
Smart-timer	10 minutes
Water resistance level	IPX3

Purpose

The NEO rejuvenation device incorporates advanced plasma technology. This technology is used to disinfect the skin and combat skin diseases by eliminating surface infections. Plasma rejuvenation works by regenerating the skin: the heat generated by the plasma stimulates collagen restoration. Plasma also restructures lipid bilayers and forms temporary channels, which facilitates the penetration of medicinal and cosmetic products into the skin.

Device Components

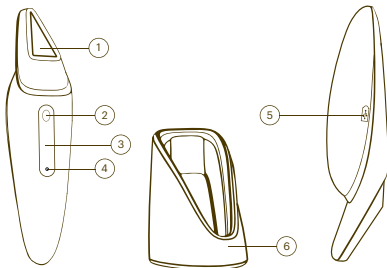1. Working Surface
2. On/ Off/ Intensity Level Switch Button
3. Intensity Level Indicator
4. Charging Indicator
5. Charging Port
6. Base

Package Contents

Device x1, Base x1, Charging Cable x1, Storage Pouch x1, Cleaning Cloth x1, User Manual x1.

Operating Procedure

1. Turn on the device by pressing and holding the on/ off button for 2 seconds. A sound will indicate that the device is on.
2. The display will light up, showing the intensity level indicators. The device offers three intensity levels: STRONG, MEDIUM, and SOFT.
3. Select the desired intensity level by pressing the level switch button repeatedly. The SOFT level will activate first, and its icon will light up on the display. Pressing the button again will switch to the MEDIUM level, and a subsequent press will activate the STRONG level.
4. The device features a smart timer; it will automatically turn off after 10 minutes of operation, accompanied by an audible signal.
5. To turn off the device, press and hold the on/ off/ intensity level switch button.

Usage Method

1. Wash your face with your preferred cleanser to remove impurities and excess sebum.
2. Dry your face thoroughly before using the massager. Any remaining water on your skin may cause significant discomfort during device application.
3. Turn on the device and select the desired intensity level. It's recommended to start with the lowest setting to allow your skin to adapt to the apparatus.
4. Slowly move the massager across your skin along the massage lines, moving upwards from the chin to the hairline. Pay special attention to problem areas. During use, the device will emit a plasma generation sound, and you may also notice an ozone smell — this is normal.
5. After completing the procedure, turn off the device, cleanse your face, and apply your skincare product. Gently massage your face to improve cosmetic absorption.
6. Clean the working surface of the device with a damp (alcohol-soaked) wipe.

[Important]: Dust and sebum residue easily adhere to the ceramic part of the device's working surface when broken down by plasma during the procedure. This may lead to yellowing of the ceramic surface but does not affect the device's performance.

Cleaning & Maintenance

- Before cleaning, make sure the device is turned off.
- Do not immerse the device in water or wash it under running water.
- Please clean the device with a soft cloth or cotton pad.
- It is permissible to clean the working surface and the body with a damp cloth.
- Do not use alkaline cleaning agents, alcohol-based essential oils, or solvents to clean the device.

Storage

- Keep the device in a dry and dust-free place.
- Keep the device away from high temperatures.
- Keep the device away from aggressive liquids.
- Keep the device out of reach of children.

Charging

- When the battery level is low, the charging indicator will flash rapidly and the device will emit short beeps before turning off.
- Connect the Type-C charging cable to the device's port, and the USB end of the cable to the adapter.
- The charging indicator flashes during charging.
- The device will not turn on while charging.
- When fully charged, the charging indicator will remain solid.

Precautions

Do not use the device:

- on children;
 - on individuals suffering from allergies;
 - on individuals who have undergone facial/ eye surgery.
-
- Do not use the device if it has been immersed in water, is malfunctioning, or if it has been damaged or broken.
 - Do not use the device in the eye area (along the upper and lower orbital bone), on the temples, on the thyroid gland, or on the chest and genital areas.
 - Do not disassemble the device, as this may lead to damage, malfunctions, electric shock, fire, or injury.

Warranty

The NEO Plasma Device by WAU comes with a 1-year warranty from the date of purchase. If technical defects are found within the warranty period, the device will be replaced or repaired.

This warranty does not cover malfunctions caused by improper use or transportation, or devices that have undergone self-repair.

NEO

PLASMA DEVICE

Model: **BPW-002**

Materials: **ABS plastic, PC plastic, silicone**

Dimensions (HxWxD): **189x45x50 mm**

Weight: **183 g**

Voltage: **3.7 V**

Power: **5 W**

Made in: **PRC**

Authorized manufacturer representative: I.E. Umantsev
A.N., Leningradsky Prospekt 37, Office 14-14, Moscow
125167, Russia.

info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com