wau

AI

EYE MASSAGER

USER GUIDE

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
Technical specifications

Parameter	Value
Model Number	B-016
Materials	ABS Plastic + Zinc Alloy
Dimensions (H×W×D)	Device: 155×26×28 mm Charging Station: 62×55×59 mm
Weight	Device: 94 g Charging Station: 67 g
Built-in Power Source	Lithium Battery
Power Supply Voltage	3.7 V
Operating Voltage	11 V
Charging Time	2 hours
Battery Capacity	450 mAh
Maximum Continuous Use Time	Up to 40 minutes
Modes	EMS + Red Light / RF + Red Light

Purpose

Thank you for purchasing our AI Eye Skin Device. The device utilizes RF and EMS technologies in combination with red light therapy (620 nm – 630 nm) to help reduce wrinkles around the eyes, dark circles, and under-eye puffiness. The device helps to tighten and soothe the skin around eyes, brighten the complexion, and improve the absorption of skincare products, enhancing their overall effectiveness.

Device Components

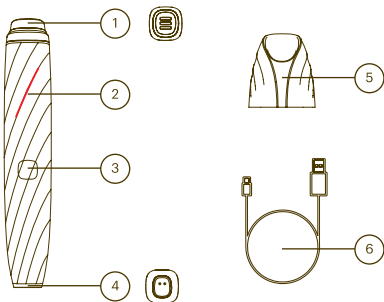- 1. Operating Surface
- 2. Mode and Battery Indicator
- 3. On/off/Mode Switch

- 4. Charging Contact Point
- 5. Charging Station
- 6. USB Type-C Charging Cable

Package Contents

Device ×1
Storage Pouch ×1
Charging Station ×1
Charging Cable ×1
User Manual ×1

Operating Procedure

1. To turn on the device, press and hold the power button.
2. The device has two operating modes: RF + Red Light Mode and EMS + Red Light Mode. Press the power button to switch between modes.

Red Indicator Light –
RF + Red Light Mode

Blue Indicator Light –
EMS + Red Light Mode

3. To turn off the device, press and hold the power button.
4. The device is equipped with an automatic shut-off feature and will power off after 10 minutes of continuous use.

Usage Instructions

Please read the instructions carefully before using the device.

[Note]: Users with reduced sensitivity to heat should use the RF function with caution, as it generates a warming effect.

RF + Red Light Mode

Step 1: Press and hold the power button to turn on the device, then select RF + Red Light Mode (indicated by a red light).

Step 2: After applying a conductive massage gel, place the working surface of the device against the skin and begin massaging. Move slowly along the massage lines: from the inner to the outer corner of the eye, following the orbital bone, and then outward in a fan-like motion from the outer corner, as shown in the diagram. Switch to the other eye every 2 minutes.

[Note]: This mode is recommended for use in the evening, no more than 4 times per week.

The radiofrequency waves generated by the device heat the dermis. Please keep the device moving during use and avoid lingering in one area to prevent burns. The device emits an audible signal every 10 seconds to remind you to change position.

The radiofrequency waves generated by the device penetrate into the subcutaneous tissues, where natural resistance produces thermal energy. The device operates on the principle that collagen fibers in the dermis contract at temperatures between 55–70°C. This tightens the skin and stimulates collagen regeneration. The immediate contraction of collagen fibers promotes skin renewal, increases dermal thickness and density, and helps reduce the appearance of wrinkles.

Using red light at a specific wavelength enhances blood circulation, has anti-inflammatory properties, reduces swelling and discomfort, brightens pigmentation, and is suitable for all skin types. Red light penetrates deeper than other wavelengths, stimulating fibroblasts – cells responsible for collagen production – which promotes skin rejuvenation and combats signs of aging.

EMS + Red Light Mode

Step 1: Press and hold the power button to turn on the device, then select EMS + Red Light Mode (indicated by a blue light).

Step 2: After applying a conductive massage gel, place the working surface of the device against the skin and begin massaging. Move slowly along the massage lines: either along the orbital bone or the brow line, from the inner to the outer corner of the eye, as shown in the diagram. Switch to the other eye every 2 minutes.

[Note]: This mode is recommended for use in the morning, no more than 2–3 times per week. The EMS microcurrents stimulate facial muscles and motor nerves, improving blood circulation and activating cellular function. This helps reduce puffiness, tighten the skin, and increase its firmness and elasticity.

Using red light at a specific wavelength enhances blood circulation, has anti-inflammatory properties, reduces swelling and discomfort, brightens pigmentation, and is suitable for all skin types. Red light penetrates deeper than other wavelengths and stimulates fibroblasts — cells that synthesize collagen — promoting skin rejuvenation and helping to fight signs of aging.

Precautions

- 1) Do not expose the charging cable to water. Do not submerge the device in water.
- 2) Do not charge the device in locations where it may get wet, such as in the bathroom, shower, or near a pool.
- 3) Do not throw the device into fire.
- 4) Use only the original charging cable provided with the device to charge it.
- 5) To avoid injury, stop using the device if any accessories or the charging cable are damaged.
- 6) If the charging cable is damaged, replace it with an original cable to Fully charge the device at least once every two months to extend battery life.
- 7) This device is not intended for use by children under 12 years old, or by individuals with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they are supervised.
- 8) Do not use the device if it is visibly damaged. Do not attempt to disassemble or repair the device yourself. Doing so may result in serious eye or skin injuries, electric shock, or even death.
- 9) Use the device only as intended and in accordance with this manual.
- 10) It is recommended to charge, operate, and store the device at temperatures between 5°C and 35°C.
- 11) Do not use the device if you have undergone cosmetic procedures in the past two weeks (e. g., chemical or mechanical peels, laser resurfacing, fractional laser treatment, mesotherapy, dermal fillers, or any other facial contouring procedures).

Recommendations

- 1) Regular cleaning of the device ensures optimal performance and extends its lifespan.
- 2) Do not use the device while sleeping or in a drowsy state.

Contraindications

- 1) If you have any medical conditions, consult a dermatologist before using the device.
- 2) Do not use the device if you are sensitive to light.
- 3) Do not use the device if you have cracks or open wounds on the skin, skin diseases or severe acne, sunburn, skin infections, etc.
- 4) Do not use the device on the upper or lower eyelids, eyeballs, or mucous membranes.
- 5) Do not use the device on moles, warts, ingrown hairs, or irritated skin.
- 6) Do not share the device with others, as this may result in cross-contamination and the spread of harmful infections.

Cleaning

Following these instructions is essential to ensure your device functions properly.

Failure to do so may reduce its effectiveness or cause it to stop working. Ensure the device is turned off before cleaning.

- a. Carefully clean the device housing and working surface after each use to extend its lifespan.
- b. Do not rinse the device or charging station under running water or submerge them in water.
- c. Use a dry cloth to clean the device.
- d. Do not use metal sponges, abrasive cleaners, or harsh liquids such as alcohol or acetone.

[Note]: Wipe all components with a clean, dry cloth.

[Tip]: You may use wet wipes to remove skincare product residue from the working surface.

Storage

- a. Store the device in a dry, dust-free place.
- b. Keep the device away from high temperatures.
- c. Keep the device away from corrosive or aggressive liquids.
- d. Store the device out of reach of children.

Charging

1) Charging via computer: turn on your computer, connect the charging cable to a USB port on the computer and to the Type-C port on the charging station. Place the device onto the charging station to begin charging.

2) Charging with an adapter: connect the USB charging cable to a USB adapter and the Type-C port on the charging station. Place the device onto the station to begin charging. Adapter specifications: 5V 1A

3) Indicator light: While charging, the indicator light will blink. Once fully charged, the indicator will stay on continuously.

[Note]: Before charging, make sure that the contact point and charging cable are dry and free of any moisture or skincare products. Do not charge the device in locations where it may fall into water. Please do not use the device while it is charging. The device should not remain on the charging station for more than 24 hours.

Disposal

Do not dispose of the device with household waste. Instead, please take it to a designated recycling center to help protect the environment.

Warranty

The AI Eye Massager by WAU is covered by a 1-year warranty. During the warranty period, the device is eligible for replacement or repair in the event of any technical defects. This warranty does not cover malfunctions caused by improper use, transportation damage, or unauthorized repairs performed by the user.

NEO

EYE MASSAGER

Model: **B-016**
Material: **ABS Plastic, Zinc Alloy**
Dimensions (H×W×D): **155 × 26 × 28 mm**
Weight: **94 g**
Voltage: **3.7 V**
Power: **1.5 W**
Made in: **PRC**

Authorized manufacturer representative: I. E. Umantsev A. N.
Leningradsky Prospekt 37, Office 14-14, Moscow 125167, Russia
info@wauglobal.com

Important Notice: To ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com

RoHS