wau

ATAMA

HAIR GROWTH & SCALP
MASSAGE DEVICE

USER GUIDE

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RoHS

Table of Contents

Technical specifications.....	22
Purpose	22
Device Components	23
Operating Procedure	24
Usage Method	25
Precautions	26
Cleaning & Maintenance.....	27
Storage.....	27
Charging	27
Disposal.....	28
Warranty.....	28

Technical specifications

Parameter	Value
Model Number	BP216
Materials	ABS Plastic + Silicone + Steel SUS303
Dimensions (H×W×D):	234x87,4x47,7mm
Weight:	213,5 g
Built-in Power Source	AC100-240V 50/60HZ
Power Supply Voltage	5V1A
Operating Voltage	3,7 V
Battery Capacity	850 mAh
Charging Time	3 hours
Operating time from full charge	2 hours
Smart-timer	Up to 10 minutes
Water Resistance Level	IPX6
Irridation Intensity	Red light — 35,659 mW/cm ² Blue light — 80,265 mW/cm ²

Purpose


The ATAMA massage device is a perfect choice for accelerating hair growth, maintaining hair health, and improving scalp condition at home.

Each user has access to several modes for individual adjustment:

EMS — stimulates the scalp with ions and microcurrents, helps to clear clogged pores and activates «sleeping» hair follicles. As a result, the scalp absorbs nutrients more easily, and the hair roots are strengthened.
LED — red/ blue light stimulates the growth of hair follicles, improves their immunity, accelerates hair growth, makes them thicker, fights hair loss, and also increases their elasticity and adds shine.

VR — microvibration massage stimulates the activation of dormant hair follicles and the growth of healthy hair.

Device Components

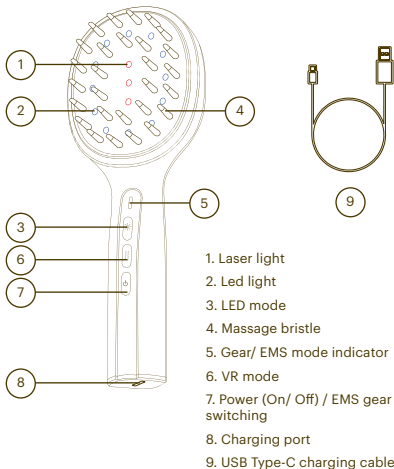Package Contents

Device ×1
Storage Pouch ×1
Charging Cable ×1
User Guide ×1

Operating Procedure

1. The massage device's working surface is protected by two protective films to prevent damage. Remove both films before using the device.

2. To turn on the device, press and hold the power button for 1,2 seconds. Upon activation, the device will emit one beep and enter EMS Mode (Level 1), indicated by one illuminated LED. A short press after activation switches to EMS Mode Level 2 (two LEDs will light up). Another short press switches to EMS Mode Level 3 (three LEDs will light up). A subsequent press turns off EMS Mode, and all three LEDs will blink simultaneously. To power off the device, press and hold the button again for 1,2 seconds. The device will emit two beeps and shut down.

3. Vibration Mode: First short press: activates continuous vibration. Second short press: activates pulsating vibration. Third short press: turns off vibration mode. Each mode switch is accompanied by a short beep.

4. LED Light Therapy Mode: First short press: activates red light therapy (630 nm) + laser (650 nm). Second short press: switches to blue light (466 nm). Third short press: turns off LED mode. Each mode switch is accompanied by a short beep.

After 10 minutes of continuous use, the device will automatically power off without entering standby mode. This means there is no remaining time buffer for switching functions.

When the battery is low, the device will emit five beeps. Just before shutting down, the battery indicator will flash ten times.

EMS, LED, and Vibration modes can be used in any combination.

Usage Method

1. Turn on the massager.

2. Select your desired combination of modes.

Note: If you are using EMS mode, ensure that your hair and scalp are slightly damp to improve conductivity. Vibration and light therapy modes are effective even on dry hair.

3. To massage the crown of the head, move the device in circular, clockwise motions around the crown area. Spend 1-2 minutes massaging each area.

4. To massage other parts of the head, move the device in the direction of hair growth – from the forehead toward the back of the head – as if combing the hair. You may also use spiral or zigzag movements.

5. Repeat all movements in each area 3 to 5 times.

6. After completing the massage, turn off the device.

Precautions

1. Please read this user guide carefully before using the device.
2. Use with caution in the following cases:
 - During menstruation or lactation;
 - In the presence of telangiectasia (spider veins);
 - When wearing temporary dental fillings or undergoing dental treatment;
 - If you have a metal allergy;
 - In cases of fever or hypertension;
 - If you have undergone plastic surgery, or are experiencing skin inflammation or itching.
3. Do not use near areas with metal implants (e. g., eyeballs, inner ear, bones, or joints).
4. Do not use on areas affected by venous or arterial disorders, thrombophlebitis, or blood clotting issues.
5. Avoid applying strong physical pressure to the scalp.
6. Use the device only as instructed in the user manual.
7. Do not disassemble, service, or modify the device yourself.
8. Store the device in a dry and clean place.
9. Do not store the device in areas with extremely high or low temperatures or near open flames.
10. Unplug the device after charging is complete.


Cleaning & Maintenance

- a. Make sure the device is turned off and not connected to a power source before cleaning.
- b. After each use, wipe the device with a damp cloth and store it in a cool, dry place.
- c. This device is intended for personal use only; do not share it with others.
- d. Keep this user guide for future reference.

Storage

- a. Store the device in a dry and dust-free place.
- b. Keep the device away from high temperatures.
- c. Keep the device away from aggressive liquids.
- d. Store the device out of reach of children.

Charging

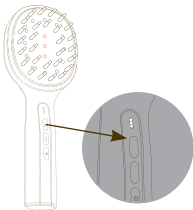1) Charging via computer connection:
Turn on the computer, connect the USB charging cable to the USB port of the computer and the Type-C port of the device.

2) Charging via adapter:
Connect the USB charging cable to the USB port of the adapter and the Type-C port of the device.

3) Indication:
During charging, 3 indicator lights will light up one by one, when 100% charge is reached, all 3 lights will be on continuously.

Power supply specification: 5V 1A 110-220V (Failure to comply with this condition will result in damage to the device and its breakdown).

Disposal

Do not dispose the device with household waste. Instead, please take it to a designated recycling center to help protect the environment.

Warranty

The ATAMA hair growth & scalp massage device by WAU is covered by a 1-year warranty. During the warranty period, the device is eligible for replacement or repair in the event of any technical defects. This warranty does not apply to malfunctions caused by improper use, transportation damage, or unauthorized repairs performed by the user.

ATAMA

HAIR GROWTH & SCALP MASSAGE DEVICE

Model: **BP216**

Material: **ABS Plastic + Silicone + Steel SUS303**

Dimensions (H×W×D): **234×87,4×47,7 mm**

Weight: **213 g**

Voltage: **5 V**

Power: **3 W**

Made in: **PRC**

Authorized manufacturer representative: WAU Global Limited
63A College road, Northcote, Auckland, 0627, New Zealand
e-mail: info@wauglobal.com

Important Notice: To ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com

RoHS