wau

LED MASK 3.0

WIRELESS PHOTOTHERAPY MASK

USER GUIDE

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Table of Contents

Technical Specifications	39
Purpose	41
Device components	41
Operating Procedures.....	43
Usage Method	44
Precautions	48
Cleaning	50
Storage.....	51
Charging	51
Disposal	52
Warranty	52
FAQ	53

Technical Specifications

Parameter	Value
Model	FM-23
Materials	ABS plastic, silicone
Dimensions	255x580x121 mm
Weight	265 g
Power source	Li-ion battery
Battery capacity	850 mAh
Operating voltage	3,7 V
Supply voltage	5 V
Smart-timer	15 minutes
Charging time	Mask — 4 hours Remote — 45 minutes
Max. continuous use time (Level 3)	Red light — 64 minutes Blue light — 60 minutes Yellow light — 55 minutes Green light — 55 minutes White light — 57 minutes NIR (Near Infrared light) — 93 minutes DIR (Deep Infrared light) — 215 minutes
Irradiation intensity	Level 1 Red — 10,5 mW/cm ² Blue — 12,4 mW/cm ² Yellow — 10,4 mW/cm ² Green — 8,4 mW/cm ² White — 13,3 mW/cm ² NIR — 7,76 mW/cm ² DIR — 3,14 mW/cm ²


Technical Specifications

Parameter	Value
Irradiation intensity	Level 2 Red — 14,4 mW/cm ² Blue — 18,5 mW/cm ² Yellow — 13,9 mW/cm ² Green — 11,2 mW/cm ² White — 16,9 mW/cm ² NIR — 10,49 mW/cm ² DIR — 4,24 mW/cm ²
	Level 3 Red — 17,4 mW/cm ² Blue — 20,2 mW/cm ² Yellow — 17,6 mW/cm ² Green — 14,6 mW/cm ² White — 19,1 mW/cm ² NIR — 17,1 mW/cm ² DIR — 8,3 mW/cm ²
Modes	Automatic and manual control
Number of LEDs	66 lamps

Purpose

The LED MASK 3.0 is a wireless 3D mask with an anatomical fit, designed for comprehensive home care. The mask features 7 key light wavelengths and 3 intensity levels to address a wide range of skin concerns. Light waves of different lengths penetrate to various depths of the dermis, triggering natural regeneration processes: Red light fights wrinkles; Blue light reduces inflammation; Green light helps with hyperpigmentation; Yellow light improves microcirculation and skin drainage; White light promotes the absorption of cosmetics; Near Infrared (NIR) heals skin damage; Deep Infrared (DIR) provides an «inner glow» effect. Every parameter of the LED MASK 3.0 (light, time, automatic mode switching) can be controlled via the remote control or buttons on the mask, creating the ideal personalized skincare routine.


Device components

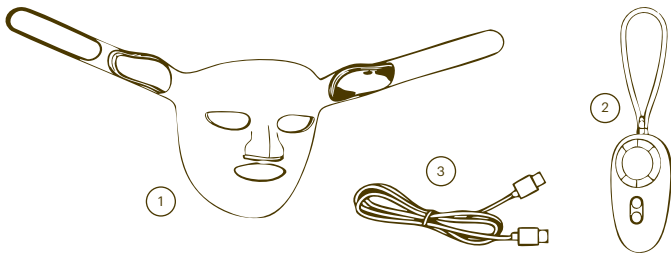1. LED MASK 3.0
2. Remote control
3. USB Type-C charging cable


Symbols on the control panel

- 4. Color selection buttons
- 5. Intensity and time adjustment buttons
- 6. Power button/Auto and Manual mode switch
- 7. Start/Pause button
- 8. Session time in minutes
- 9. Remote control battery level
- 10. C (Color) Light type:
R (red) — Red light
G (green) — Green light
B (blue) — Blue light
RG (red + green) — Yellow light
RGB (red + green + blue) — White light
IR_ — Near Infrared
IR^ — Deep Infrared
- 11. L (Level) Intensity level: 1/2/3
- 12. Operating mode: A (auto) — Automatic control
M (manual) — Manual control

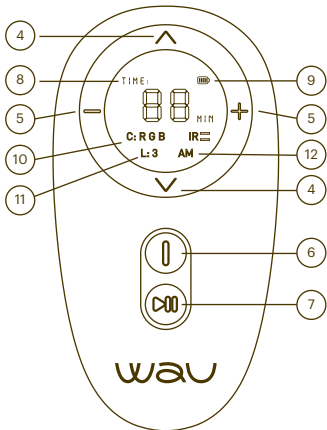Markings on the mask controller

- 13. Power button/Color selection
- 14. Intensity adjustment/
Language selection

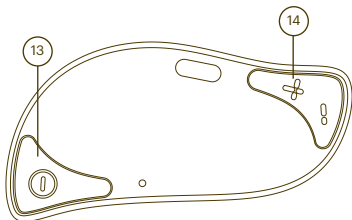Package contents

Mask x1, remote control x1, USB Type-C charging cable x1, user manual x1, storage bag x1.

Operating Procedures

LED MASK 3.0 features a dual control system:

1. Via the controller on the mask body.
2. Via remote control.

Control via the built-in controller

When using the built-in controller, the mask can only operate in automatic mode.

1. Press and hold the button for 2 seconds to turn the mask on or off. A long beep will sound to confirm.
2. Short press the button to switch between light colors.
3. Press and hold the button for 5 seconds to enable Bluetooth and connect to the remote control. Once the mask is connected to the remote, the indicator on the mask controller will turn blue.

1. Short press the button to increase the light intensity level.
2. Press and hold the button to select the voice prompt language: English or Russian.

1. Short press the button to decrease the light intensity level.
2. Press and hold the button to turn the voice prompts on or off.

Control via the remote control

To control the device via the remote, first turn on the mask using the built-in controller, then turn on the remote and wait for the mask to connect to the remote control. Upon successful connection, a voice prompt will sound: «Bluetooth is connected/Bluetooth подключен».

1. Press and hold the button to turn the mask on.
2. Short press the button to select the control mode: Automatic (one light color operates throughout the entire session, selected by the user) or Manual (the user sets a sequence of colors, defining the duration and intensity for each).
3. Press and hold the button to turn the mask off.

1. Short press to confirm the selected settings: intensity level and time.
2. Double press after configuring all settings in Manual mode to start the mask operation.
3. Short press during the session to pause the treatment. Short press again to resume the session. The screen flashes while paused.
4. In Manual mode, press and hold after entering the mode to reset the preset settings (the last configured settings).

Short press the button to navigate to the desired light: the «down» button switches to the next light type; the «up» button returns you to the previous light settings.

1. In Manual mode, short press the increase button to increase the intensity level and session time; press the decrease button to reduce them.
2. In Automatic mode, short press the increase button to increase the light intensity; press the decrease button to decrease the intensity level.

Usage Method

Sensitivity test

Important: You must carry out a sensitivity test before your first treatment.

1. Place the device on the inside of your forearm.
2. Turn on the mask and run at least a 10-minute phototherapy session.
3. If you feel heat or discomfort, stop using the device immediately.
4. After completing the test, wait at least 6 hours. If the skin on the inside of your arm is red and/or itchy, the LED MASK 3.0 is not suitable for you. Contact the organization authorized to handle consumer claims.

After completing the sensitivity test

1. Cleanse your face before starting the session: remove makeup, do not apply skincare products.
2. Turn on the mask using the controller.
3. Place the mask on your face so that it fits snugly.
4. Connect and fasten the Velcro straps at the back so the mask sits comfortably.

Automatic control mode

5. The automatic mode can be controlled either directly via the mask or via the remote control.

To control directly via the mask

6. Turn it on by pressing and holding the Power button on the controller.

7. Short press the same button (Power/Mode switch) to select the desired light type.

The following light types are available in LED MASK 3.0:

Red Light type (630 nm)
Remote control symbol — R

Stimulates collagen synthesis and effectively fights wrinkles.

Blue Light type (415 nm)
Remote control symbol — B

Provides an antibacterial effect and reduces inflammation.

Yellow Light type (630 nm + 520 nm)
Remote control symbol — RG

Improves microcirculation and provides skin lymphatic drainage.

Green Light type (520 nm)
Remote control symbol — G

Evens out skin tone and helps with hyperpigmentation.

White Light type (630 nm + 520 nm + 415 nm)
Remote control symbol — RGB

Accelerates metabolism and promotes deep absorption of cosmetics.

Near Infrared (NIR) (830 nm)
Remote control symbol — IR₋

Accelerates regeneration and heals micro-damage.

Deep Infrared (1070 nm)
Remote control symbol — IR⁺

Activates cellular energy, providing a pronounced «inner glow» effect.

8. Use the increase and decrease buttons to set the desired light intensity level.
9. The session will start automatically; there is no need to confirm the start.
10. A 15-minute smart timer will automatically end the program.

To control via the remote control

11. You must establish a Bluetooth connection between the mask and the remote. After turning on the mask, turn on the remote control and wait for the beep.
12. The automatic control mode is activated first on the remote. Check that the letter A is displayed in the lower corner of the screen.
13. Select your desired light.
14. Adjust the intensity level from 1 to 3.
15. The mask will start automatically, and the smart timer countdown will begin.
16. After 15 minutes, the mask will turn off automatically.

Manual control mode

17. Turn on the mask using the controller.
18. Manual mode can only be operated via the remote control; therefore, it is important to establish a Bluetooth connection. After turning on the mask, turn on the remote and wait for the beep.
19. Select the manual control mode on the remote. The letter M will light up on the screen.
20. You can independently set the desired combination of lights, intensity, and duration.
21. Red light is the first to be configured. You can start by setting the intensity level. The intensity indicator flashes on the screen. Use the +/- buttons to select the level. Confirm it by pressing the Start/Pause button.
22. Next, you can select the duration for this light. Use the +/- keys to increase or decrease the time. Confirm this setting by pressing the Start/Pause button.

23. Proceed to configure the next colors in the same way (Blue → Yellow → Green → White → Near Infrared → Deep Infrared).
24. If you wish to skip a color, set its operating time to zero or skip it by moving to the next light using the V/Λ navigation buttons.
25. Using the Λ button, you can return to the previous light to adjust its brightness and duration.
26. Double-click the start/pause button to begin the session.
27. The LED MASK 3.0 features a memory function in manual mode. Pressing and holding the Start/Pause button after entering manual settings on the remote will reset the preset light combination.
28. After the phototherapy session, it is recommended to apply skincare and moisturizing products.

Frequency and duration of use

The mask is equipped with a smart timer; the session will end automatically after 15 minutes. It is not recommended to use the device longer than the specified time. If you are using the phototherapy mask for the first time, we recommend following the schedule below:

Week 1: Every 3 days for 5–15 minutes. Start with a 5-minute session and gradually increase the duration.

Week 2: Every 3 days for 15 minutes.

Week 3 and onwards: Every 2 days for 15 minutes.

To achieve a visible effect, it is recommended to use the LED MASK 3.0 for at least 4 weeks.

Precautions

LED MASK 3.0 is a device for home use. It must only be used in accordance with the method specified in this user manual. Use for any other purposes is prohibited, as it may result in product damage or physical harm.

This device is not intended for use by individuals (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the safe use of the device by a person responsible for their safety.

Contraindications

- Do not use the mask if you suffer from lupus erythematosus, photosensitive eczema, or albinism. Using phototherapy under these conditions may cause a severe skin reaction.

- It is prohibited to use the mask if you suffer from any photosensitivity disorder (if you are sensitive to light).

- Do not use the mask if you are taking medications that cause photosensitivity. These may include certain antibiotics, chemotherapy drugs, and diuretics. Consult your physician before using the device.

- Do not use the mask if you have open wounds, or avoid the affected area during the phototherapy session.

- Do not use the mask if you have had a sunburn within the last 2 weeks.

- Do not use the mask if you have Fitzpatrick skin type I, V, or VI.

- Do not use the mask if you have a history of light-triggered epilepsy seizures.

- Do not use in cases of severe acne.

- It is prohibited for infants, people with thyroid diseases, or those taking thyroid medications to use the device.

- Do not use the mask during charging or maintenance. Do not touch the mask, charging cable, or USB adapter with wet hands while charging.

- Do not use the mask if it is damaged. Contact the organization authorized by the manufacturer to handle consumer claims.

Warnings

- A sensitivity test must be carried out before using the mask (see «Instructions for use»).
- Turn off and remove the mask if it feels hot or if your skin feels uncomfortably warm.
- Use the mask only on cleansed skin: remove makeup and do not apply skincare products before the procedure. Some ingredients may cause photosensitivity or reduce the treatment effect.
- Keep the mask away from water and avoid getting water on the device.
- Do not place the mask on hot surfaces.
- Do not attempt to modify or repair the device yourself. This will void the warranty and may result in electric shock. The device contains no user-serviceable parts.
- Do not use the mask if the device itself, the controller, the charging plug, or the cable is damaged.
- Before connecting, check if the mains voltage matches the voltage specified on the device.
- Do not subject the equipment to strong impacts to avoid damaging the lamps.
- Do not look directly at the light source while the device is operating to protect your eyes.
- Do not use the mask in the following environments: near MRI, diathermy, electrocoagulation, wireless power transfer, 5G, or other systems with electromagnetic transmitters.
- Avoid using this equipment adjacent to or stacked with other equipment, as it could result in improper operation. If such use is necessary, both devices should be observed to verify they are functioning normally.
- Use original converters and cables. Using accessories other than those specified or provided by the manufacturer may result in increased electromagnetic emissions or decreased immunity of this equipment, leading to improper operation.
- Do not use the mask closer than 30 cm to operating radio devices: Wi-Fi routers, mobile phones, baby monitors, and antennas. Strong radio signals may cause malfunctions or reduce the effectiveness of the procedure.

Miscellaneous

- LED MASK 3.0 is a personal hygiene item; it is not recommended to share the device with other persons.
- If the mask is damaged, stop using it immediately and contact the organization authorized to handle consumer claims to avoid any danger. In the event of self-repair, disassembly, or incorrect use, the authorized organization shall be released from all warranty obligations.
- Immediately stop using the device if you experience discomfort or any adverse reaction (skin redness, itching, allergy) as a result of using the mask.

Cleaning

- Clean the inner surface of the mask after each use.
 - Turn off the mask before cleaning.
 - Do not use high-temperature or aggressive disinfectants.
 - Do not use benzene, thinners, butyl acetate, or other flammable substances for cleaning.
 - Do not spray water directly onto the surface of the mask.
1. Surface cleaning:
 - Wipe both sides of the mask with a soft, lint-free cloth dampened with a mild water-based cleaning solution.
 - Place the mask on a clean, dry surface until it is completely air-dried.
 2. Visual inspection:
 - Inspect the device surface for any visible dirt. If contamination remains, proceed to the next step.
 3. Deep cleaning (for stubborn stains):
 - Carefully wipe the silicone surface twice with a 70%–80% ethanol solution (3 minutes each time), then wipe thoroughly dry with a damp or dry cloth.
 - Allow the mask to air-dry completely before use.

Storage

1. Away from flammable and explosive materials (sprays, oils, gas, etc.).
2. At room temperature in a dry place.
3. Avoid storing in dusty areas, in direct sunlight, or in high-humidity conditions.

Charging

To operate the LED MASK 3.0, both the controller on the mask and the remote control must be charged.

1. Charging the mask

- Low battery ($\leq 3,5V$): After turning on, the indicator on the controller turns red and flashes slowly.
- Fully discharged ($\leq 3,3V$): The indicator flashes rapidly 5 times, then the mask turns off automatically. It cannot be turned on until recharged; attempts to turn it on will result in the indicator flashing 5 times again.
- Ensure the mask is turned off.
- Connect one end of the Type-C cable to the controller's charging port.
- Connect the other end to a 5V 1A power adapter; the red indicator on the controller will light up.
- During charging, the indicator will stay lit continuously.
- When fully charged, the indicator will turn blue.
- Full charging time is approximately 4 hours.

2. Charging the remote control

- Ensure the remote is turned off.
- Connect the Type-C cable to the remote's charging port and a 5V 1A power adapter; a battery icon will appear on the screen.
- Charging progress is shown by the filling bars of the icon.

 0-25% заряда — One bar flashes.

 25-50% заряда — One bar stays lit, second flashes.

 50-75% заряда — Two bars stay lit, third flashes.

 75-100% заряда — Three bars stay lit, fourth flashes.

 100% заряда — All four bars stay lit while connected. The display turns off once disconnected.

Important

- Use only the original USB Type-C cable included in the package.
 - We recommend an adapter with at least 0.5 A current for faster charging.
 - The device features overvoltage protection.
-
- Normal parameters: Voltage 5V, current 35 mA to 1 A.
 - Protection trigger: If a 9V or 12V adapter is used, the device will block the charging function.
 - Risk of damage: If the charging voltage exceeds 16V, the device may be permanently damaged.

Disposal

At the end of the device's service life, its components (especially the lithium-ion battery and electronic parts) may harm the environment if disposed of with regular household waste.

1. Do not dispose of the device in a regular waste bin along with food waste.
2. Dispose of the device in accordance with local legal regulations.
3. Take the device to specialized electronics and battery collection points for further sorting and recycling.

Warranty

The WAU brand LED MASK 3.0 Wireless Phototherapy Mask is covered by a 1-year warranty. During the warranty period, the product is subject to replacement or repair if technical defects are discovered. This warranty does not cover malfunctions caused by improper operation or transportation, nor does it apply to products after unauthorized repair.

FAQ

Q: Is this mask waterproof?

A: No, this mask is not waterproof. Washing the mask is prohibited. For more detailed maintenance instructions, please refer to the «Cleaning» section.

Q: How long can the mask operate on a full charge?

A: The operating time for each light at intensity level 3 (the most energy-consuming level) is specified in the «Technical Specifications» section.

Q: My skin feels dry after using the mask. How can I fix this?

A: Since the light wavelengths affect the dermis and stimulate cell metabolism, increasing the skin's absorption capacity, it requires more hydration and nourishment. Therefore, after using the LED MASK 3.0, it is recommended to replenish the skin's moisture levels by applying sufficient hydrating and nutrient-rich products.

Q: Does the mask heat up during use?

A: When using Red, Blue, or Green light, the mask temperature does not increase. When using combined lights (Yellow, White), the temperature may rise to 38 ± 1 °C.

Q: Will I see the effect immediately after the first use?

A: The effect after one or two sessions will be less noticeable because the skin's metabolism follows a 28-day cycle. Therefore, visible results will typically appear after 1 month of consistent use.

Q: What is the recommended session duration?

A: Phototherapy should not exceed 15 minutes per session; otherwise, it may harm the skin. Please monitor your usage time.

Q: Can I bend or fold the mask?

A: No, as this may lead to the failure of its internal circuitry.

LED MASK 3.0

WIRELESS PHOTOTHERAPY MASK

Model: **FM-23**

Material: **ABS plastic, silicone**

Dimensions (H×W×D): **255x580x121 mm**

Weight: **265 g**

Voltage: **3,7 V**

Power: **3 W**

Origin: **PRC**

Authorized manufacturer representative: WAU Global Limited
63A College road, Northcote, Auckland, 0627, New Zealand
e-mail: info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com