wau

MAMI

EMS & PHOTOTHERAPY MASSAGER

USER GUIDE

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RoHS

Table of Contents

Technical Specifications	26
Purpose	26
Device Components	27
Operating Procedure	28
Usage Method	29
Cleaning & Maintenance	31
Storage	31
Charging	32
Precautions	33
Disposal	34
Warranty	34




Technical Specifications

Parameter	Value
Model	ST-H313
Nominal Operating	5-9 V
Nominal Power Supply Voltage	3,7 V
Power Source	Lithium-ion Battery
Battery Capacity	900 mAh
Material	ABS Plastic, Zinc Alloy
Dimensions	Device: 147,1×83×44,4 mm Base: 57,1×57,1×76 mm
Weight	Device: 179 g Base: 108 g
Charging Time	3 hours
Operating Time	Red/ Blue Light Mode: 200 min Mixed Light Mode: 300 min
Modes	EMS (3 levels) + Phototherapy (3 light modes) + Vibration (optional)
Irridation Intensity	Red light — 7,7867 mW/cm ² Blue light — 32,758 mW/cm ² Mixed (Purple) light — 23,322 mW/cm ²
EMS	Level 1 — 0,08 mA Level 2 — 0,11 mA Level 3 — 0,13 mA

Purpose

The MAMI lifting massager for the face and neck by the WAU brand is a device that uses advanced microcurrent EMS therapy technology. Three levels of EMS intensity, combined with low-frequency vibration and red, blue, and mixed LED light therapy, improve the skin's absorption ability and help tighten it. The result of using MAMI is healthy, elastic, and radiant skin.

Device Components

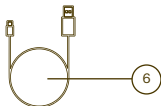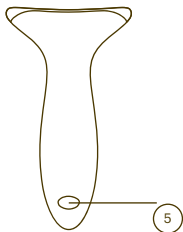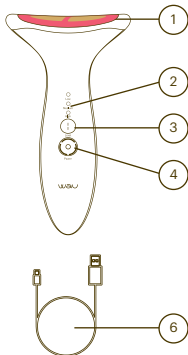- 1. Massage Surface
- 2. EMS Intensity level Indicator/
Charging Indicator
- 3. LEVEL Button — Adjusts
EMS intensity/ Turns vibration
on or off


- 4. POWER Button — Power on/
off/ Mode switch button/ Mode
indicator
- 5. Charging Port
- 6. USB Type-C Charging Cable
- 7. Base

Package Contents

Massager ×1
Base ×1
USB Charging Cable ×1
User Guide ×1
Storage Pouch ×1
Cleaning Cloth ×1

Operating Procedure

1. This device features 3 phototherapy light modes and 3 levels of EMS intensity.
2. To turn on the device, press and hold the POWER button. The POWER button also allows you to select the phototherapy light mode.
3. After the massager is turned on, use the LEVEL button to select the EMS intensity level from low to high (short press) or to switch the vibration function on/ off (long press).
4. To turn off the device, press and hold the POWER button.

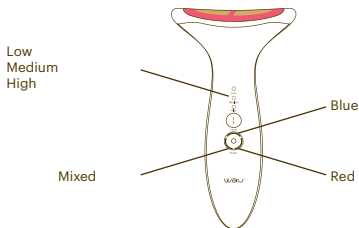
Smart Timer: The device will automatically turn off after 10 minutes of operation.

Usage Method

Step 1:

Press and hold the POWER button to turn on the massager. Use short presses of the POWER button to select the phototherapy light mode. The device will automatically start with the blue light. Press the POWER button again to switch to red light. Press the POWER button once more to activate the mixed light mode.

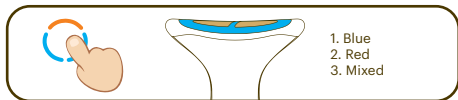Step 2:

Use short presses of the LEVEL button to select the EMS intensity level. The intensity levels range from low to medium, and then to high.

Low level: 1 indicator lit (Low)

Medium level: 2 indicators lit (Low + Medium)

High level: 3 indicators lit (Low + Medium + High)


Press and hold the LEVEL button to turn the vibration function on or off.

Step 3:

Apply a conductive massage gel to the skin. Place the massage surface of the device on the skin of your face.

Step 4:

Massage your face along the lines shown in the illustration, moving from the bottom (chin) upwards (towards the hairline on the forehead). Repeat the massage along the indicated lines at least three times. It is not recommended to use the device for more than 10 minutes at a time.

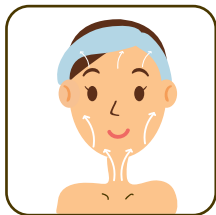Effects and Features of Each Mode:

1. EMS acts on muscle fibers, causing them to actively contract and relax. EMS therapy strengthens the muscular framework, increases tone, and helps tighten the skin by improving muscle function. Intensive stimulation enhances blood circulation, allowing tissues to receive more oxygen and nutrients. The overall effect of EMS use is visibly stronger muscles, improved contours, and increased skin firmness.

2. Phototherapy

Blue Light: Fights acne, brightens the skin, tightens pores, soothes sensitive skin, and regulates sebum production.

Red Light: Stimulates blood and lymphatic circulation and activates mitochondria within cells. The specific wavelength of red light converts light energy into cellular energy, neutralizing free radicals in the cells. Red light also triggers fibroblast production, which generates collagen, promoting improved blood circulation, wrinkle reduction, and deep skin hydration.

Mixed Light: Combines the benefits of both red and blue light.

3. Vibration

Activating vibration massages and firms the skin, reducing the appearance of fine lines and wrinkles.



Note: Use the MAMI massager up to 4–5 times per week for best results.

Cleaning & Maintenance

Please carefully follow the instructions below to ensure proper maintenance of the device. Failure to adhere to these recommendations may reduce the massager's effectiveness and performance.

Before cleaning, make sure the device is turned off

- Use only dry cloths to clean the device.
- Do not use metal scrubbers, abrasive cleaning agents, or corrosive liquids such as alcohol or acetone.
- Thoroughly clean the device and check its functionality after each cleaning.
- Do not rinse the device under running water or immerse it in water.
- Ensure all parts of the device are completely dry before connecting it to the charger.

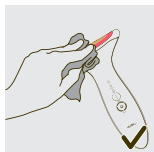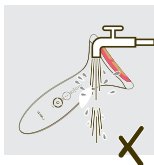Storage

- Store the device in a dry, dust-free environment.
- Keep the device away from high temperatures.
- Avoid exposure to aggressive liquids.
- Store the device out of reach of children.

Charging

1. Charging via computer using USB cable:

Turn on your computer and connect the USB cable to the computer's USB port. Connect the other end of the cable with the Type-C connector to the device to start charging.

2. Charging with an adapter:

Plug the USB cable into a power adapter and then connect the Type-C end of the cable to the device's charging port.

Adapter specifications: 5V = 1A.

3. During charging, all three microcurrent intensity indicators will flash. When fully charged, the indicators will remain steadily lit.

Attention: before charging, ensure the charging port on the device is dry and free from water or cosmetic residues.

Please do not charge the massager in locations where it may come into contact with water. Do not use the device while it is charging.

Charging time should not exceed 24 hours.

To extend the battery life, fully charge the massager at least once every two months.

Please use only the cable provided with the massager.

Precautions

Please read the user manual carefully before using the massager and keep this instruction for future reference.

Danger

- Keep the charging port and USB cable dry. Never immerse them in water.
- Do not charge the device in a humid environment (e.g., near a filled bathtub, running shower, or a pool).
- Do not throw the device into fire.

Caution

- Use only the USB cable supplied with this massager.
- To avoid injury, do not use the massager or USB cable if they are broken or damaged.
- If the USB cable is damaged, replace it only with an original cable to avoid hazards.
- Fully charge the massager at least once every two months to prolong battery life.
- This device is not intended for use by children under 12 years old, people with limited physical, sensory, or mental capabilities, or those without experience or knowledge of such devices, unless they are supervised or instructed on proper use.

Use with Care

- Use the device only as intended.
- Charge, use, and store the massager at temperatures between 5 °C and 40 °C (35 °F to 95 °F).
- Do not use the device if your skin is in the healing phase after cosmetic procedures such as chemical peels, laser resurfacing, fractional laser treatments, mesotherapy, or dermal filler injections within the last two weeks.

General Recommendations

- Regular cleaning of the device ensures effective performance and extends its service life.
- Do not use the massager while sleeping or in a drowsy state.

Additional Warnings

- Do not apply the device directly to the ear. This may cause temporary hearing loss.
- Do not use this device if you have ear health issues (such as otitis externa, ear infections, or other ear diseases).
- Do not use the massager on cracked skin, open wounds, or problematic or irritated skin (such as skin affected by severe acne, sunburn, infection, etc.).
- Do not use the device if you have a pacemaker, implants, or metal prostheses.
- Do not use the device on the sensitive skin around the eyes, as this skin is thin and delicate.
- Do not use the device on mucous membranes (e.g., inside the nose, mouth, or genital area), as this may cause local irritation.
- Do not use this device if you have thyroid gland disorders or are taking medication for thyroid treatment.

Disposal

Do not dispose of this device with regular household waste. Please dispose the massager in accordance with local regulations or take it to a designated recycling center.

Warranty

Thank you for purchasing the MAMI EMS & Phototherapy massager by WAU. This device is covered by a 1-year warranty from the date of purchase. During the warranty period, products found to have technical defects are subject to repair or replacement.

This warranty does not cover malfunctions caused by improper use, transportation damage, or devices that have been repaired independently.

MAMI

EMS & PHOTOTHERAPY MASSAGER

Model: **ST-H313**

Material: **ABS plastic, zinc alloy, stainless steel**

Dimensions (HxWxD): **147x83x44 mm**

Weight: **179 g**

Voltage: **3,7 V**

Power: **5 W**

Made in: **PRC**

Authorized manufacturer representative: WAU Global Limited
63A College road, Northcote, Auckland, 0627, New Zealand
info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com

RoHS