wau

MIRA

FACE AND NECK LIFTING DEVICE

USER GUIDE

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RoHS

Table of Contents

Technical Specifications	24
Purpose	25
Device Components	25
Operating Procedure	26
Usage Method	28
Cleaning	29
Storage	29
Charging	30
Precautions	31
Disposal	31
Warranty	32


Technical Specifications

Parameter	Value
Model	MR-991
Materials	PC plastic, ABS plastic, zinc alloy, silicone
Dimensions	179,5x72x48,4 mm (device) 81x69,8x64,8 mm (base)
Weight	164 g (device) 115 g (base)
Rated Power	10 W
Operating Voltage	3,7 V
Power Source	Built-in lithium polymer battery
Battery Capacity	1350 mAh
Modes	Custom combinations
Technologies	Microcurrents (3 intensity levels) Phototherapy 3 LED lights: Red — 630 nm, Blue — 460 nm, Green — 520 nm Wave Vibration (7000 rpm) Heating (42±3 °C) Cooling (15±3 °C)
Smart-Timer	10 minutes *Heating mode: 5 minutes *Cooling mode: 3 minutes
Operating Time	35 minutes
Charging Time	3 hours
Waterproof Level	IPX4

Purpose

The MIRA multifunctional facial device combines 5 technologies, including 3 levels of microcurrents, 3-color LED therapy, wave vibration, and heating and cooling modes for comprehensive skin rejuvenation, allowing you to freely combine modes to create a personalized skincare routine featuring microcurrents for cellular energy and contour lifting, LED therapy for photo-rejuvenation and targeting imperfections, vibration for lymphatic drainage and reducing puffiness, heating for deep nourishment and muscle relaxation, and cooling for pore tightening and instant skin soothing.

Device Components

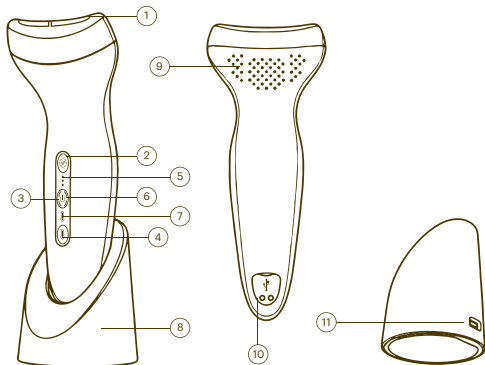1. Treatment surface
2. Microcurrents/Vibration On/Off
3. Power On/Off/LED therapy switch button
4. Cooling/Heating On/Off
5. Microcurrent level indicator
6. Phototherapy mode indicator
7. Cooling/Heating mode indicator
8. Charging base
9. Ventilation holes
10. Type-C port for direct charging and pins for contactless charging
11. Type-C port for base charging

Package Contents

Massager x1, charging base x1, treatment surface cover x1, USB Type-C charging cable x1, storage pouch x1, user manual x1.

Operating Procedure

1. **Power On:** Press and hold the middle button (Power On/Off) to turn on the device. An activation signal will sound.
2. **Standby Mode:** After powering on, the massager enters standby mode, and the LED therapy indicators will pulse gently. If no mode is selected within 2 minutes of entering standby, the device will shut down automatically.
3. **LED Therapy:** Briefly press the middle button to cycle through the LED light modes. The first press activates red light, the second – blue, the third – green, and the fourth press turns the phototherapy function off.
4. **Visual Feedback:** With each LED light switch, a light-emitting diode around the center button will light up.

5. **Microcurrents:** Briefly press the top button (Microcurrents On/Off) to activate the microcurrent mode. Short presses allow you to set the intensity level: first press for Level 1 (low), second for Level 2 (medium), and third for Level 3 (high). An additional press will turn the microcurrent mode off.
6. **Intensity Display:** Each intensity level is shown on the display: one LED for Level 1, two LEDs for Level 2, and three LEDs for Level 3.
7. **Wave Vibration:** Press and hold the top button (Microcurrents/Vibration) to activate wave vibration. Press and hold again to turn it off.
8. **Cooling Mode:** Briefly press the bottom button (Cooling/Heating) to activate cooling mode. It operates for 3 minutes and shuts off automatically to prevent overheating and maintain the lowest possible temperature. A blue thermometer icon will appear on the display. Please note: if the device is turned off during cooling mode or if the smart timer expires, the cooling fan will continue to run for 10 seconds before shutting down automatically.
9. **Heating Mode:** Press the temperature button again to switch from cooling to heating. A red thermometer icon will appear on the display. Heating mode operates for up to 5 minutes before automatic shutdown.
10. **Manual Extension:** If you wish to continue using cooling or heating after the automatic shutdown, simply restart the desired mode.
11. **General Smart Timer:** The device is equipped with a 10-minute master smart timer; it will shut down automatically after 10 minutes of total operation.
12. **Power Off:** To manually turn off the device, press and hold the middle Power On/Off button.

Usage Method

1. Preparation: Apply a conductive gel or a water-based gel skincare product to your cleansed face and neck. This ensures optimal conductivity for the microcurrents to achieve maximum results.
2. Power On: Turn on the device and select your desired modes.

3. Neck Massage: Start with the neck area. Move in vertical strokes downward from the earlobe toward the collarbones to "open" the lymphatic pathways.

4. Facial Massage: Massage the face by moving along the massage lines from the center outward toward the periphery. Direction of movement: diagonally upward to counteract gravity and lift the muscles. Perform the procedure moving from bottom (chin) to top (hairline).

5. Forehead: To smooth horizontal forehead wrinkles, move the device in vertical upward strokes.

6. Preparation Stage (up to 2 minutes): It is recommended to start with Heating and Vibration modes. This combination prepares the skin, making it highly receptive to active cosmetic ingredients and stimulation. The goal is to relax muscles, expand pores for cleansing, and boost blood circulation and lymphatic drainage.

7. Main Treatment Stage: This stage stimulates cells and muscles while improving blood flow to promote collagen and elastin production, muscle toning (lifting), and targeting specific skin concerns (e.g., acne or pigmentation). Activate Microcurrent mode at a comfortable intensity level and turn on LED Therapy, starting with the Red light.

8. Targeted Treatment (total stage duration — 6 minutes): Following the initial rejuvenation, you can switch from Red light to Blue light (to treat acne and oiliness) or Green light (to soothe the skin and address rosacea, redness, or couperose).

9. Finishing Stage (up to 2 minutes): Complete the massage using the Cooling function. This stage locks in the results. The cold will soothe the skin, reduce residual redness, and "seal" the active ingredients of your skincare products within the skin.

Cleaning

- Turn off the device and disconnect the charging cable before cleaning.
- Wipe the body and the treatment surface with a soft, damp cloth.
- Do not use organic solvents such as alcohol or gasoline.
- Do not immerse the device completely in water.

Storage

- If the device is not used for an extended period, store it fully charged and recharge the battery every 3 months.
- Store in a cool, dry place, avoiding direct sunlight and high-temperature environments.

Charging

- Use only the original or recommended power adapter for this device.
- Place the device in a well-ventilated area during charging; do not cover it with any objects.
- If you notice abnormal heating, smoke, or other unusual phenomena during charging, disconnect the power immediately.
- Do not use the device while it is charging.

1. Charging via Base: Place the massager into the charging base. Ensure that the device contacts are firmly aligned with the base pins. Connect the USB end of the charging cable to a USB adapter and the Type-C end to the charging station.

2. Direct Charging: Connect the USB end of the charging cable to a USB adapter and the Type-C end directly to the device after opening the silicone cover.

3. Indicator: During charging, the 4 indicators on the Power button will light up sequentially in a circular pattern.

4. Full Charge: When 100% charged, all 4 indicators will remain lit.

Precautions

- This device is intended for personal home use only and is not a medical device.

• The following individuals should use the device only under the guidance of a physician:

- Pregnant and lactating women.
- People with sensitive skin or skin diseases.
- People with implanted electronic medical devices (e.g., pacemakers).
- Diabetics or people with circulatory disorders.

• Do not use on the following areas:

- Damaged skin (wounds, inflammation, eczema).
- Mucous membranes (eyes, mouth).

• When using Cooling mode, do not keep the device on the same skin area for too long to prevent frostbite.

• When using Heating mode, stop use immediately if you feel excessive heat to prevent burns.

• Children must use the device only under adult supervision.

Disposal

Do not dispose of this product with household waste. Please take it to a specialized recycling center to help protect the environment.

Warranty

The MIRA Face and Neck Lifting Device by WAU is covered by a 1-year warranty from the date of purchase. If technical defects are found within the warranty period, the device will be replaced or repaired. This warranty does not cover malfunctions caused by improper use or transportation, or devices that have undergone self-repair.

MIRA

FACE AND NECK LIFTING DEVICE

Model: **MR-991**

Materials: **PC plastic, ABS plastic, silicone, zinc alloy**

Dimensions: **179,5x72x48,4 mm (device),**

81x69,8x64,8 mm (base)

Weight: **164 g (device), 115 g (base)**

Voltage: **3,7 V**

Power: **10 W**

Origin: **PRC**

Authorized manufacturer representative: PERFECT SQUARE
LIMITED 54/2 Fairclough Rd, Beach Haven, Auckland,
Auckland, 0626, New Zealand
e-mail: info@wauglobal.com

Important Notice: to ensure continuous product improvement,
the manufacturer reserves the right to modify design and
technical specifications. The latest version of this manual
is available at: wauglobal.com