wau

MIGA BLACK

ORAL CARE IRRIGATOR

USER GUIDE

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Table of Contents

Technical Specifications	25
Purpose	25
Device Components	26
Operating Procedure	27
Usage Method	28
Cleaning & Maintenance	29
Storage	30
Charging	30
Troubleshooting and Solutions	31
Precautions	32
Disposal	33
Warranty	33


Technical Specifications

Parameter	Value
Model Number	RM-MINIO18
Nozzle Diameter	0,6 mm
Pulsation Rate	1600 pulses per minute
Water Pressure	60-140 PSI
Water Tank Capacity	140 ml
Battery Capacity	1100 mAh
Noise Level	≤ 72 dB
Input Voltage	5 V = 1 A
Rated Power	8 W
Charging Port Type	Type-C
Charging Time	4 hours
Operating Time	90 minutes
Water Resistance Rating	IPX7
Modes	3


Purpose

The MIGA Black oral care irrigator by WAU is a portable device designed for effective dental hygiene. Equipped with three operating modes and four interchangeable nozzles, it efficiently removes food particles, dental plaque, and cleans interdental spaces and gum pockets. The device features a compact design, making it convenient for travel and ideal for maintaining oral hygiene with orthodontic appliances. The irrigator offers long battery life and complies with FDA standards.

Device Components

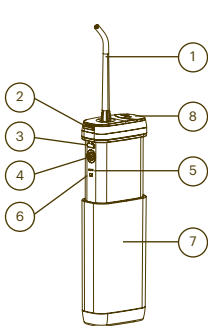- 1. Multifunctional Nozzle
- 2. Nozzle Release Button
- 3. Charging Port
- 4. Power Button
- 5. Mode Indicator
- 6. Battery Indicator
- 7. Water Tank

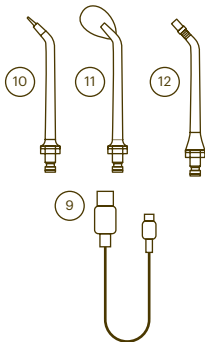- 8. Nozzle Storage and Water Filling Compartment
- 9. Charging Cable
- 10. Periodontal Nozzle
- 11. Tongue Cleaning Nozzle
- 12. Orthodontic Nozzle

Package Contents

Irrigator ×1
Storage Pouch ×1
Nozzles ×4
Charging Cable ×1
User Manual ×1

Operating Procedure


1. Insert the nozzle until it clicks into place.

2. Open the water tank and fill it with clean water.

3. Place the nozzle in your mouth and briefly press the power button to turn the device on.

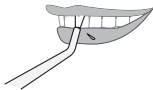4. To switch between modes, press and hold the power button. After use, turn off the device with a short press. The irrigator will remember the last used mode and automatically start in that mode the next time it is turned on.

5. To remove the nozzle, press and hold the nozzle release button and gently pull it out. Open the storage compartment and place the nozzle inside.

6. Remember to clean the device regularly: detach the water tank and rinse it thoroughly.

[Note]: We recommend starting with the gentle mode, especially during the initial period of use. Slight bleeding may occur during the first few uses — this is normal and should subside within a week. If bleeding persists for more than two weeks, discontinue use and consult your dentist.

Usage Method

a. Use the standard multifunctional nozzle. Press the power button and select the desired mode. Position the nozzle so that the water jet is perpendicular to the area being treated.

[Note]: Sensitive gums may bleed easily. Do not direct the water jet directly at the gums!

b. Begin cleaning from the back teeth and move toward the front teeth. Clean the entire dental arch by pausing briefly at each tooth until the gum line and interdental spaces are thoroughly cleaned. If you started on the outer side of the teeth, repeat the procedure on the inner side. Perform this routine for both the upper and lower teeth.

[Note]: For people with sensitive teeth, it is recommended to use the gentle mode.

c. To clean gum pockets, use the periodontal nozzle. Press the power button and select the appropriate mode. Position the nozzle at an angle so that the tip is directed at the gum line and into the periodontal pocket. Clean according to the instructions described above.

d. For tongue cleaning, use the spoon-shaped nozzle. Attach the nozzle, press the power button, and select the required mode. Move the nozzle along the tongue with sweeping motions in one direction, from the root to the tip of the tongue. Tongue cleaning should last about 30 seconds, which is approximately 5 passes.

e. The orthodontic nozzle is designed for cleaning braces. In addition to the water jet opening, the orthodontic nozzle features a stiff-bristle brush at the tip, which effectively cleans the brackets by removing plaque and food debris.


Indicators and Cleaning Modes

Mode	How to Switch	Indicator Color	Recommended For	Water Pressure
Gentle	Short press of the power button	Green	Low pressure; ideal for sensitive teeth and first-time users	40-60 PSI
Pulse	Long press of the power button	Blue	Alternating pulsating pressure for deep cleaning and gum massage	60-110 PSI
Powerful	Long press of the power button	White	High pressure; effective daily cleaning for most users	120-140 PSI

Cleaning & Maintenance

- After use, turn off the irrigator and remove the nozzle.
- Clean and dry all components thoroughly, and store the device in a dry place.
- This device is rated IPX7 waterproof, allowing it to be rinsed directly under running water.

Storage

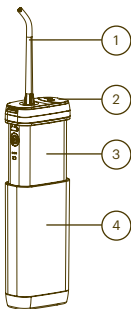1. Nozzle

Press the nozzle release button and remove the nozzle. Rinse it with water and dry thoroughly with a cloth.

It is recommended to replace the nozzle every 6 months.

If the tip of the nozzle becomes deformed, replace it immediately.

2. Water Filling Compartment

Clean the compartment regularly to prevent contamination and leakage.

3. Device Body

Wipe off dirt with a cloth dampened in mild detergent, then dry with a clean cloth. If the device is heavily soiled, you may rinse it with water. Do not use solvents, gasoline, alcohol, or similar substances to clean the device. Do not wash with hot water above 50°C. Always dry the device thoroughly after cleaning.

4. Removable Water Tank

Open and rinse the removable water tank regularly. Clean the inside with a brush or suitable cleaning tool, then wipe dry with a cloth. If the device will not be used for an extended period, do not leave water inside the tank.

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Troubleshooting and Solutions

Issue	Possible Cause	Solution
Device does not turn on	Newly purchased device or not used for more than 3 months	Charge continuously for more than 4 hours
Low water pressure	a. Nozzle is clogged with debris. b. Nozzle is deformed	Replace the nozzle
Runs only a few minutes after charging	a. Insufficient charging time, low battery. b. Battery has reached end of life	a. Charge continuously for more than 4 hours b. Contact customer service
Battery indicator does not light up during charging	a. Charging cable is not properly connected b. Power circuit is damaged	a. Check the cable connection b. Contact customer service

Precautions

1. Do not use the device without first reading the operating and safety instructions to avoid injury caused by improper use.
2. Do not disassemble the device or replace its components, as this may cause damage and void the warranty.
3. Do not use hot water (above 50°C / 122°F) in the water tank for oral cleaning to prevent burns.
4. If the power cable or plug is damaged, replace it before use.
5. Keep the device away from open flames and heat sources. Do not charge, use, or store the device in high-temperature environments.
6. Use the device only for its intended purpose: oral hygiene.
7. Do not apply excessive pressure with the nozzle on your teeth or gums to avoid injury.
8. Do not immerse the device in water for prolonged periods to prevent malfunction, fire, electric shock, or other hazards.
9. Do not damage the power cable: avoid stretching, bending, or twisting it. Do not place heavy objects on the cable or trap it between items. Keep both the device and cable away from heat sources.
10. Do not store the device in hot or humid environments, and avoid direct sunlight.
11. Do not charge the device in hot or damp places, such as bathrooms or toilets.
12. If the device will not be used for a long time, unplug the power adapter and charging cable to prevent fire or other risks.
13. This product is not intended to be repaired or refurbished. If service is required, please return it to the retailer.
14. Do not operate the device without a nozzle attached.
15. This device is not intended for use by children. Use only under adult supervision. Keep out of reach of children. Do not use it as a toy.
16. If you have gum disease or have undergone oral surgery within the past two months, consult your dentist before using the device.
17. Do not add mouthwash or any other liquids to the water tank. This may cause corrosion and the tank to stick to the device body, making it difficult to open and possibly damaging the unit.

Disposal

Do not dispose of this product with regular household waste. Instead, take it to a designated recycling center to help protect the environment.

[Battery Handling]:

The built-in battery contains substances that may harm the environment. When disposing of the irrigator, be sure to remove the battery and take it to a certified battery recycling facility.

Do not disassemble the battery itself—only remove it prior to disposal. Ensure the battery is fully discharged before removal.

[Battery Removal Instructions]:

1. Disconnect the irrigator from the power source.
2. Turn off the device.
3. Open the device casing and remove all internal components.
4. Cut the positive (+) and negative (–) battery terminals.
5. Carefully remove the battery.

Warranty

The MIGA Black oral care irrigator by WAU is covered by a 1-year warranty. During the warranty period, if any technical defects are discovered, the product is subject to replacement or repair.

This warranty does not cover malfunctions resulting from improper use or transportation, or products that have been repaired independently.

MIGA BLACK

ORAL CARE IRRIGATOR

Model: **RM-MINI018**
Material: **ABS Plastic**
Dimensions (H×W×D): **124 × 62 × 42 mm**
Weight: **183 g**
Voltage: **5 V**
Power: **8 W**
Made in: **PRC**

Authorized manufacturer representative: PERFECT SQUARE LIMITED
54/2 Fairclough Rd, Beach Haven, Auckland, Auckland, 0626,
New Zealand
info@wauglobal.com

Attention! Due to continuous product improvement, the manufacturer reserves the right to make changes to the design, components, and technical specifications.
The latest version of this manual can be found at wauglobal.com

RoHS