wau

MIO2

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## Table of Contents

---

Technical Specifications .....	26
Purpose .....	26
Device Components .....	27
Operating Procedure .....	28
Usage Method .....	28
Cleaning & Maintenance .....	30
Storage.....	30
Charging .....	31
Precautions .....	32
Disposal.....	33
Warranty.....	33

## Technical Specifications


Parameter	Value
Model	ST-F818
Materials	ABS plastic, PC plastic, zinc alloy
Dimensions (H×W×D)	Massager: 136×54×22 mm Charging base: 85×78×52 mm
Weight	Massager: 80 g Charging base: 70 g
Power Source	Lithium-ion battery
Battery Capacity	500 mAh
Operating Voltage	3,7 V
Power	5 W
Charging Time	2,5 hours
Operating Time (full charge)	50 minutes
Smart Timer	10 minutes
Modes	EMS/ Heating/ EMS + Heating + Vibration
EMS	≤0,50 mA

## Purpose

The MIO2 microcurrent device with EMS function uses EMS technology, heating, and vibration to maintain youthful, healthy, and fresh-looking skin. The device helps reduce puffiness, enhances the absorption of skincare products, and improves facial contour. Its specially curved shape is perfect for Gua Sha massage, while the pointed ends can be used for acupressure massage and muscle relaxation.

## Device Components

---

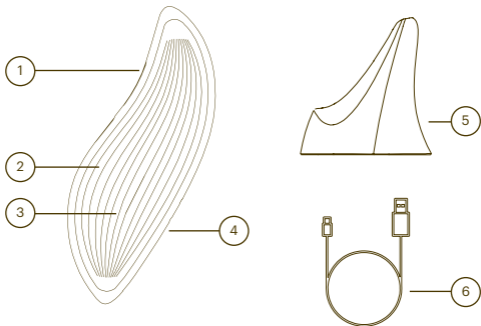1. Working surface
2. Power/ mode/ charging indicator
3. Power on/ off/ mode switch button
4. Magnetic charging port
5. Charging base
6. Charging cable

### Package Contents

Device ×1, charging base ×1, charging cable ×1, storage pouch ×1, cleaning cloth ×1, user guide ×1.

## Operating Procedure

---

1. Press and hold the power/ mode button to turn on the device.
2. The massager has 3 operating modes. After turning it on, the EMS mode starts first. The indicator lights up blue.
3. Press the power/ mode button to switch modes. The next mode is heating. During this mode, the indicator glows red.
4. Pressing the mode button again will activate the “EMS + heating + vibration” mode. The indicator lights up white.
5. To turn off the device, press and hold the power/ mode button.
6. The device features a smart timer. After 10 minutes of operation, it will automatically shut off. Note that when switching modes, the smart timer does not reset but continues counting down until the end of the session, regardless of the current mode.

## Usage Method

---

1. Turn on the device and select the desired mode.
2. Apply a conductive massage gel to the skin. Place the massaging surface of the device against the facial skin.
3. Perform the facial massage, following the massage lines.
4. Move the device upward from the chin toward the hairline on the forehead.
5. Repeat the procedure along all massage lines at least three times. It is not recommended to use the device for more than 10 minutes.

[Note]: During massage in the third mode, it is recommended to close your eyes

and relax your face. Then place the massager on areas with wrinkles or areas prone to wrinkle formation, such as crow's feet and nasolabial folds — for example, the area between the eyebrows; right/ left temples; right/ left nasolabial folds.

6. Clean the massage surface with a tissue.

### *Effects and Features of Each Mode:*

1. EMS acts on muscle fibers, causing them to actively contract and relax. Intensive muscle stimulation increases their tone and elasticity, which helps lift the skin and improve facial contours. Enhanced blood circulation during EMS therapy helps tissues receive more oxygen and nutrients. The overall effect of EMS use is strengthening the muscular framework, improving skin tone, and achieving a more pronounced lifting effect.

2. In Heating mode, gentle warmth (up to 45°C) stimulates blood flow, accelerates skin metabolism, soothes the skin, and enhances the absorption of active ingredients such as Vitamin C and minerals. This mode is recommended for use with skincare products.

3. The combination of all three technologies in the EMS + Heating + Vibration mode stimulates circulation, calms the skin, activates skin cells, reduces puffiness, boosts collagen production, and increases skin elasticity. The effects of this mode may include smoothing of fine lines and wrinkles, lifting of facial and neck contours, and the creation of a more defined V-shaped facial profile.

## Cleaning & Maintenance

---

- Please clean the device after each use, as residue from cosmetic products may damage it.
- Ensure the device is turned off before cleaning.
- Use only dry cloths to clean the device.
- Do not use metal scrubbers, abrasive cleaning agents, or corrosive liquids such as alcohol, acetone, or similar substances to clean the massager.
- Clean the device thoroughly and check its functionality after each cleaning.
- Do not rinse the device under running water or immerse it in water.
- Make sure all parts of the device are completely dry before connecting it to the charger.

## Storage

---

- Store the device in a dry, dust-free environment.
- Keep the device away from high temperatures.
- Keep the device away from harsh or corrosive liquids.
- Store the device out of reach of children.

## Charging

---


1. Charging via computer using USB cable: Turn on the computer, connect the USB cable to the computer's USB port, and connect the other end with the Type-C connector to the charging base.

2. Charging using an adapter: Connect the USB cable to a power outlet via an adapter. Connect the other end with the Type-C connector to the charging base.

Adapter specifications: 5 V, 1 A

3. During charging, the indicator light will blink. Once fully charged, the indicator will stay lit continuously.

- After a prolonged period of inactivity, it may take approximately 15 minutes for the device to reactivate and start charging again.
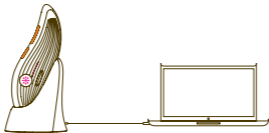
## Precautions

---

Please read the user manual carefully before using the massager and keep this instruction for future reference.

### Dangerous

- Never immerse the device in water.
- Do not charge the device in a damp environment (e. g., near a filled bathtub, running shower, or pool).
- Do not throw the device into fire.

### Warning

- Do not use the device if the massager or any of its components are damaged.
- Fully charge the massager at least once every 2 months to prolong battery life.
- This device is not intended for use by children under 12 years old or individuals with limited physical, sensory, or mental capabilities, or those without experience or awareness of such massagers, unless supervised or instructed personally.
- Do not insert metal objects into the device's charging port to avoid short circuits.

### Caution

- Use the device only as intended.
- Charge, use, and store the massager at temperatures between 5°C and 40°C (35°F to 95°F).
- Do not use the device if you are taking steroid medications.
- Do not use the device if your skin is healing from cosmetic procedures such as chemical peels, laser resurfacing,

fractional laser rejuvenation, mesotherapy, or dermal fillers within the last two weeks.

— Consult a dermatologist if you have any concerns regarding the use of this device considering your skin condition or medical contraindications.

— If you experience discomfort during use, stop immediately and consult a specialist.

— For hygiene reasons, do not share the device with others.

### **General Recommendations**

— Regular cleaning ensures effective operation and extends the device's lifespan.

— Do not use the massager while sleeping or drowsy.

### **Disposal**

---

Do not dispose of this device with regular household waste. Please dispose of the massager in accordance with local regulations or take it to a designated recycling center.

### **Warranty**

---

The MIO2 microcurrent massager with EMS function by the WAU brand is covered by a 1-year warranty from the date of purchase.

If any technical defects are discovered during the warranty period, the product is eligible for repair or replacement.

This warranty does not cover damage caused by misuse, improper transportation, or unauthorized repairs.

# MIO2

---

## MICROCURRENT MASSAGER WITH EMS FUNCTION

---

Model: **ST-F818**

Materials: **ABS plastic, PC plastic, zinc alloy**

Dimensions (H×W×D): **136×54×22 mm**

Weight: **150 g**

Voltage: **3,7 V**

Power: **5 W**

Made in: **PRC**

Authorized manufacturer representative:

WAU Global Limited 63A College road, Northcote,

Auckland, 0627, New Zealand

[info@wauglobal.com](mailto:info@wauglobal.com)

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: [wauglobal.com](http://wauglobal.com)

RoHS