wau

MIGA

ORAL CARE IRRIGATOR

USER GUIDE

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
Technical Specifications

Parameter	Value
Model	OI-2107
Rated Supply Voltage	5 V = 4 W
Rated Operating Voltage	3.6 V = 8 W
Power Source	Lithium Battery
Battery Capacity	1050 mAh, 3.6 V
Water Resistance Rating	IPX7
Modes	Gentle/ Strong/ Pulse
Pulsation Frequency	1600 pulses per minute $\pm$ 20%
Water Tank Capacity	130 ml
Operating Time	Gentle Mode – 50 min Strong Mode – 45 min Pulse Mode – 100 min Sterilization – 220 min
Charging Time	3 hours


Purpose

The MIGA irrigator by WAU is a portable oral care device. It features 3 operating modes and a nozzle sterilization mode, an enlarged liquid reservoir, and comes with 3 attachments to effectively remove food debris, dental plaque, and clean interdental spaces and gum pockets. The device has a foldable design, making it convenient to take on trips or use for orthodontic appliance care.

Device Components

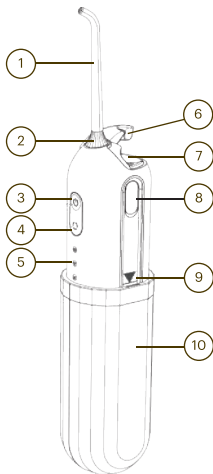- 1. Nozzle
- 2. Nozzle socket
- 3. Power button
- 4. Mode switch button
- 5. Mode indicators
- 6. Dust cover
- 7. Nozzle storage compartment

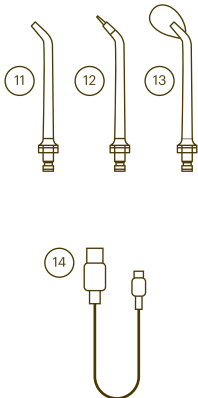- 8. Type-C charging port
- 9. Indicator
- 10. Water tank
- 11. Standard nozzle for daily tooth cleaning
- 12. Periodontal tip for cleaning gum pockets
- 13. Tongue cleaning tip
- 14. Type-C charging cable

Operating Procedure

Installing and Removing the Nozzle

[Note]

If any debris is found during the first use of the device (in the water tank or the water supply tube), it may appear that the product is damaged. Do not worry if you notice traces of clear liquid inside the device — this is normal. The liquid is distilled water used to test the irrigator's functionality. Rinse the product before use and proceed with confidence.


a. Remove the nozzle from the nozzle storage compartment.

b. Insert the nozzle into the nozzle socket — the magnetic lock will secure it automatically.

c. To remove the nozzle, gently pull it out of the socket.

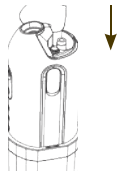d. Place the nozzle back into the storage compartment and ensure it is securely fixed.

[Caution] Do not apply pressure to the nozzle during operation.

Filling with Water

- Open the dust cover.
- While holding the irrigator vertically, fill the tank with water.
- After filling, close the dust cover securely.
- It is recommended to use clean tap water. If the water feels too cold, use room-temperature water (not exceeding 40°C). You may also use saltwater or mouthwash, but only as recommended by your dentist.

Powering On and Selecting a Mode

To switch between modes, press the Mode Switch Button:

- First indicator: Gentle Mode
- Second indicator: Strong Mode
- Third indicator: Pulse Mode
- Fourth indicator: Sterilization Mode

Select the appropriate mode.

Usade Method

- Use the standard multifunctional nozzle. Press the Power Button, select the desired mode, and position the nozzle in front of your teeth so that the water stream is perpendicular to the treatment area.

[Note]: Sensitive gums may bleed easily. Do not aim the water jet directly at the gums!

- Start cleaning from the back teeth and move forward to the front teeth.

Clean the entire dental arch, pausing briefly at each tooth until the gumline and the space between teeth are thoroughly rinsed. If you start with the outer side of the teeth, repeat the process on the inner side.

Repeat the procedure for both upper and lower rows of teeth.

[Note]: People with sensitive teeth are advised to choose Gentle Mode.

c. During use, hold the device vertically to prevent the magnetic nozzle lock from disconnecting, which may cause the water flow to stop.

d. To clean gum pockets, use the periodontal tip.

Press the Power Button and select the desired mode.

Position the nozzle at an angle, so the tip is directed at the gumline and into the gingival sulcus (gum pocket).

Perform the cleaning as described in the instructions above.

e. To clean the tongue, use the spoon-shaped tongue cleaning tip.


Attach the tip, press the Power Button, and select the appropriate mode.

Move the tip across the tongue with sweeping motions in one direction — from the root to the tip of the tongue.

Tongue cleaning should take 30 seconds, approximately 5 passes.

Sterilization Mode

a. Place the nozzle into the nozzle storage compartment and close the dust cover.

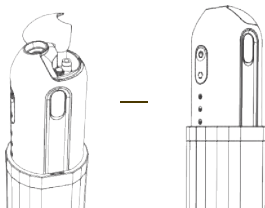b. Press the Mode Switch Button to select the fourth mode — Sterilization Mode. Then press the Power Button.

c. The sanitization process will begin, sterilizing the nozzle and the storage compartment using UV light.

Cleaning

After Use

- Turn off the device and open the dust cover.
- Turn the device upside down and close the cover.
- This device has IPX7 water resistance, so it can be rinsed directly under running water.
- Wash thoroughly, dry completely, and store in a dry place.

Storage

a. Nozzle

Rinse the nozzle with water and wipe it dry with a cloth. It is recommended to replace the nozzle every 6 months.

If the nozzle tip becomes deformed, replace it immediately.

b. Body

Wipe the body with a soft textured cloth dampened with a small amount of mild detergent.

c. Tube

The tube can be washed with or without detaching it from the device. To detach, gently pull it out. After cleaning, reattach the tube and wipe it dry. Do not bend the tube excessively, as this may cause damage.

d. Filter

The filter should be cleaned regularly, ideally every three months.


e. Water Tank

The water tank is removable for easy cleaning.

Recommended cleaning frequency: once a week.

[Recommendation]

After using the irrigator, detach and dry the water tank occasionally to maintain hygiene

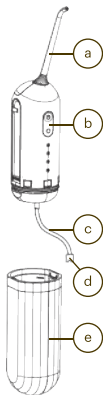[Note]: Do not use hot water above 50°C for cleaning.

[Recommendation]: Do not use thinner, alcohol, or other harsh cleaning agents, as they may cause discoloration or damage to the device.


Charging

[Note]: When the battery is low, the indicator around the button will blink red.

- Turn off the device and retract the water tank.
- Insert the Type-C charging cable into the charging port.
- Plug the adapter into an outlet and ensure that the indicator light turns on.
- While charging, the indicator around the button blinks. Once fully charged, the indicator stays solid.

A full charge takes approximately 3 hours.

With a full charge, the device can operate for around 100 minutes.

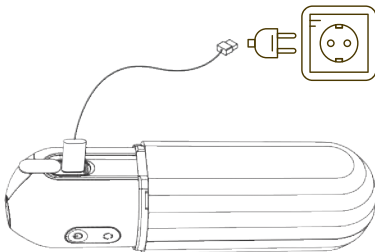[Note]: While charging, place the device on a flat surface and ensure that no water enters the water tank or charging port.

Malfunctions and Solutions

Malfunction	Cause	Solution
Does not work	New device purchased, or device has not been used for more than 3 months	a. Charge continuously for 3 hours until fully charged. b. Contact customer service.
Works only for a few minutes after charging	a. Insufficient charging time, low battery level. b. Battery life has expired.	a. Charge the device continuously for 3 hours until fully charged. b. Contact customer service.
Low water pressure	a. Nozzle clogged with dirt. b. Nozzle tip is deformed. c. Dirt in the pump.	a/b. Replace the nozzle. c. Contact customer service.
No water stream	a. Water tank is empty. b. During use, the irrigator is tilted too much. c. Water tube is improperly connected d. Dirt in the pump.	a. Fill the water tank. b. Use the device holding it vertically. c. Connect the water tube properly. d. Contact customer service.

Precautions

1. Keep the charging port dry during charging.
2. Do not connect the device to a power source with a voltage different from the specified one.
3. Do not touch the charger with wet hands.
4. Do not use the device if it is damaged or malfunctioning.
5. Do not aim the water jet directly at the tongue, ears, or other parts of the body — water pressure may cause injury.
6. If you have dental conditions, use this device under the supervision of a dentist.
7. This device can be used together with mouthwash or other products recommended by your dentist.
8. Do not disassemble the device or insert foreign objects into any openings unless you are a qualified technician.
9. Do not operate the device in harsh outdoor environments.
10. Do not fill the water tank with iodine, bleach, or water-insoluble aromatic oils, as these substances may damage the device and shorten its lifespan.
11. Remove all oral jewelry before using the device.
12. Do not use the device if there are visible wounds on the tongue or inside the mouth.
13. If prescribed antibiotics for treatment, do not use this device unless specifically advised by your doctor.
14. Do not direct the water jet into the ears to avoid injury.
15. Keep the device out of reach of children under 3 years old.
16. Store the device in a dry, clean place away from direct sunlight and heat sources.
17. Before storing the device for a long period (over 3 months), fully charge the battery.
18. Do not operate the device without water in the tank, as this may cause malfunction.

[Notes]:

1. Do not use hot water above 50°C to clean the water tank. To prolong the life of the tank, avoid deforming it due to exposure to high temperatures. It is recommended to use room temperature water below 40°C or running tap water for cleaning this device.
2. Do not add toothpaste, salt, hard mouthwash tablets, or other chemicals to the water to avoid clogging and damaging the device.
3. Before adding mouthwash or fresheners, check that they do not contain isopropyl myristate, which can damage the device and affect its normal operation.
4. Do not use thinners, alcohol, or other harsh cleaning agents, as this may cause discoloration and malfunction of the device.

Disposal

Do not dispose of the device with household waste. Instead, take it to a specialized recycling center to protect the environment.

[Battery Handling]

The built-in battery contains substances that may pollute the environment. When disposing of the irrigator, be sure to remove the battery and recycle it at an official recycling center.

Do not disassemble the battery; only remove it prior to disposal. Before removing the battery, ensure it is fully discharged.

[Battery Removal]

1. Disconnect the irrigator from the power source.
2. Turn off the irrigator.
3. Open the housing and remove all parts.
4. Cut off the positive and negative battery terminals.
5. Remove the battery.

The name and content of dangerous substances in the products

Part Name	Toxic and hazardous substances or elements					
	Cadmium (Pb)	Mercury (Hg)	Chromium (Cd)	Hexavalent chrome [Cr(VI)]	Polybrominated biphenyl (PBBs)	Polybrominated biphenyl (PBDes)
Plastic part	○	○	○	○	○	○
Hardware	○	○	○	○	○	○
Circuit Control Assemblies	×	○	○	○	○	○
Motor assembly	×	○	○	○	○	○
Battery	○	○	○	○	○	○
Cable	×	○	○	○	○	○

This table is prepared in accordance with the requirements of SJ/T 11364.

○ : indicates that the content of hazardous substances in all homogeneous materials of the component is below the maximum allowable limit specified in GB/T 26572.

×

 : indicates that the content of hazardous substances in at least one homogeneous material of the component exceeds the maximum allowable limit specified in GB/T 26572.

The service life for environmental protection according to this marking is 10 years. During this period, under normal conditions, the product containing toxic and hazardous substances or elements will not leak or suffer structural damage. Non-use of the product by the user will not cause serious environmental pollution or significant harm to health or property. Disposal of the product after the environmental protection service life should be conducted responsibly and in accordance with applicable national laws and regulations.

Warranty

The MIGA Oral Care Irrigator by WAU is covered by a 1-year warranty. During the warranty period, if technical defects are found, the product is subject to replacement or repair.

This warranty does not cover malfunctions caused by improper use, transportation damage, or products that have been repaired by unauthorized persons.

MIGA

ORAL CARE IRRIGATOR

Model: **OI-2107**
Material: **ABS plastic**
Dimensions (HxWxD): **145 x 55 x 54 mm**
Weight: **195 g**
Voltage: **5 V**
Power: **8 W**
Made in: **PRC**

Authorized manufacturer representatie: I.E. Umantsev A.N.
Leningradsky Prospekt 37, Office 14–14, Moscow, 125167, Russia
info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com

RoHS