wau

LED MASK 2.0

PHOTOTHERAPY FACE MASK

USER GUIDE

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Table of Contents

Technical Specifications	26
Purpose	27
Device Components	27
Operating Procedure	29
Usage Instructions	29
Safety Precautions	31
Cleaning	31
Storage.....	32
Charging	32
Disposal.....	33
Warranty.....	33
Frequently Asked Questions.....	33


Technical Specifications

Parameter	Value
Article Number/SKU	K9C-1KB
Materials	ABS plastic, silicone
Dimensions (H×W×D)	190x679x5 mm
Weight	390 g (inc. remote)
Power Source	Lithium-ion battery
Battery Capacity	2450 mAh
Operating Voltage	3,7 V
Input Voltage	5 V
Charging Time	3–4 hours
Maximum Continuous Usage Time	100 minutes — red/blue/green light mode 78 minutes — yellow light mode 60 minutes — white light mode
Modes	red/blue/green/yellow/white
Number of LEDs	117 lamps
Irradiation Intensity	Level 1 Red light — 6,595 mW/cm² Blue light — 17,112 mW/cm² Green light — 6,612 mW/cm² Yellow light — 14,049 mW/cm² White light — 27,542 mW/cm² Level 2 Red light — 14,058 mW/cm² Blue light — 30,635 mW/cm² Green light — 15,636 mW/cm² Yellow light — 17,857 mW/cm² White light — 32,169 mW/cm² Level 3 Red light — 16,368 mW/cm² Blue light — 47,653 mW/cm² Green light — 23,095 mW/cm² Yellow light — 28,081 mW/cm² White light — 39,089 mW/cm²

Purpose

The WAU LED MASK 2.0 is an advanced skincare device utilizing phototherapy technology. Featuring five distinct light wavelengths, it effectively addresses various skin concerns. The light therapy stimulates collagen regeneration, which enhances subcutaneous metabolism and restores skin elasticity. With regular use, this rejuvenating phototherapy mask helps maintain skin tone, combat visible signs of aging, and improve skin imperfections.

Device Components

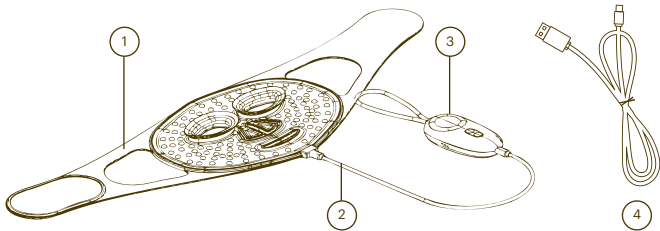
1. Phototherapy mask
2. Type-C cable for connecting the mask to the control unit
3. Control unit
4. USB-Type-C charging cable

Control Unit Indicators

1. Power On/Off Button
2. Pause Button
3. Light Intensity Adjustment Buttons
4. Mode Selection Buttons
5. Operating Time
6. Charging Indicator
7. Light Mode Indicator
8. Light Intensity Level Indicator

Package Contents

- Phototherapy mask ×1
- Control unit ×1
- Type-C cable ×1
- USB Type-C charging cable ×1
- Storage pouch ×1
- Instruction manual ×1

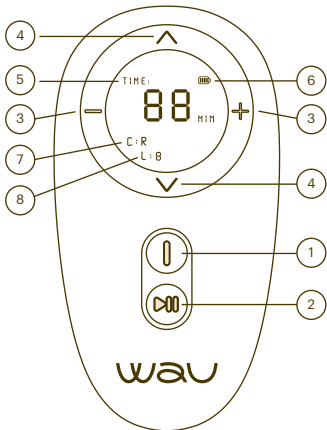Operating Instructions

1. Connect the mask to the control unit using the Type-C cable. Insert the cable into the mask's port until it clicks, then repeat the same with the control unit.
2. Power on the device by pressing and holding the On/Off button on the control unit for 2 seconds. The mask will automatically activate in red light mode at the first intensity level (control unit display indicators: C: R, L: 1).
3. Switch light modes by briefly pressing (1 second) the mode selection buttons ($\wedge/\vee$) to cycle through: Blue light (C: B); Green light (C: G); Yellow light (C: RG); White light (C: RGB).
4. Adjust intensity by briefly pressing (1 second) the light intensity buttons ($-/+$) to increase up to level 3 or decrease back to level 1 (control unit display: L: 1 (2/3)).
5. Pause operation by briefly pressing (1 second) the pause button. Press again to resume.
6. Power off the device by pressing and holding the On/Off button on the control unit for 2 seconds.

Usage Instructions

[Important] Before first use, perform a sensitivity test:

- Place the mask on the inner side of your arm.
- Turn on the mask and complete a full 10-minute phototherapy session.
- If you experience any discomfort or burning sensation on your skin, discontinue use immediately.
- Wait approximately 6 hours after testing. If the inner side of your arm becomes red and/or itchy, the LED MASK 2.0 is not suitable for you.

1. Cleanse and dry your hands and face thoroughly.
2. Secure the mask on your face by fastening the hook-and-loop straps, ensuring a comfortable fit without slipping.
3. Turn on the mask by connecting it to the control unit.
4. Select your preferred light mode and intensity level.
5. Relax and enjoy your phototherapy session.

Each light mode operates for 15 minutes before the control unit automatically switches off. When changing modes, the timer resets and begins a new 15-minute cycle.

[Important] During phototherapy, carefully monitor your total session time (when using multiple modes in one session) to ensure it does not exceed 15 minutes.

6. After completing your light therapy session, apply your preferred skincare product (serum or cream), as the skin's absorption capacity is enhanced following LED mask use.

7. For visible results, use the mask regularly 3-4 times per week for a minimum of 4 weeks.

8. After achieving your desired results, you may reduce frequency to 1-2 maintenance sessions per week.

Light Mode Specifications

Red Light (660 nm)

With the longest wavelength among the available options, red light penetrates deeper into the skin than other light modes. It targets fibroblasts, stimulating collagen production and promoting skin regeneration. As a result, red light therapy accelerates wound healing, firms and tightens the skin, and reduces visible signs of aging, including wrinkles.

Blue Light (460 nm)

Blue light has the shortest wavelength in the visible light spectrum. It targets the epidermis layer, delivering antimicrobial effects. This light is particularly effective for combating acne as it eliminates bacteria that cause breakouts.

Green Light (525 nm)

Green light is absorbed by both the epidermis and dermis, triggering a cellular response that promotes skin recovery and healing. It is particularly effective for reducing redness, soothing irritation in sensitive skin, and addressing rosacea by minimizing capillary size to improve blood flow.

Yellow Light (Red + Green = Yellow)

The yellow light in this mask version is achieved by combining red and green optical light spectra. This polychromatic light therapy delivers comprehensive effects for more visible results: yellow light stimulates blood circulation, brightens skin tone, and reduces the appearance of pigmentation.

White Light (Red + Blue + Green = White)

The blended white light spectrum penetrates the dermis and stimulates cellular metabolism, enhancing the skin's absorption capacity. Following a white light session, the skin becomes more receptive to active ingredients in skincare products. For optimal results, apply your cosmetic treatment immediately after the session to maximize its efficacy.

Safety Precautions

- The rejuvenating phototherapy mask is intended for home use only. Use the device strictly in accordance with the methods specified in this user manual. Any unauthorized use may cause product damage or bodily harm.
- If the mask is damaged, immediately discontinue use and contact authorized customer service to avoid potential hazards. Any unauthorized repairs, disassembly, or improper use will void the manufacturer's warranty.
- Do not use the mask while charging. During charging, do not touch the mask, charging cable, or USB adapter with wet hands.
- Immediately discontinue use if you experience discomfort or any adverse reaction (skin redness, itching, allergic response) from using the mask.
- The device is prohibited for use by infants, individuals with epilepsy, those with thyroid disorders, or persons taking thyroid medication.
- Do not use the mask if you have damaged, injured, or post-operative skin.

Cleaning Instructions

1. Disconnect the device from the power source before cleaning.
2. Do not spray water directly on the device during cleaning.
3. Wipe the mask surface twice with a 70-80% ethanol solution (3 minutes each time), then use a damp or dry cloth.
4. Avoid flammable cleaning agents such as benzene, solvent, or butyl acetate.
5. Do not use other cleaning products, as they may damage the mask.

Storage

1. Clean the surface of the LED MASK 2.0 after each use.
2. Store the device in a dry place at room temperature.
3. Avoid storing the device in dusty areas, under direct sunlight, or in humid environments.
4. This mask contains no user-serviceable parts.

Charging

The battery level of the LED MASK 2.0 is displayed on the control unit.

 100% charge: All 4 battery indicators are lit.

 75-100%: 3 indicators remain steadily lit.

 50-75%: 2 indicators remain steadily lit.

 25-50%: 1 indicator remains steadily lit.

 0-25%: 1 indicator flashes.

To charge the LED MASK 2.0, connect the charging cable to a USB port on your computer (or adapter) and to the Type C port on the control unit. Do not connect the charging cable directly to the mask.

The charging progress of the LED MASK 2.0 is also displayed on the control unit.

 0-25% charge: One indicator flashes.

 25-50% charge: One indicator remains steadily lit while the second flashes.

 50-75% charge: Two indicators remain steadily lit while the third flashes.

 75-100% charge: Three indicators remain steadily lit while the fourth flashes.

 100% charge: All four indicators remain steadily lit as long as the control unit is connected to the power source. The display automatically turns off when disconnected.

— If the input voltage drops below 3,2 V, the charging indicator will blink 3 times before the control unit automatically powers off.

— For charging, please use only the original USB cable included with your device. Do not connect the mask directly to the charging USB cable.

— Using an adapter with voltage below 5 V will prevent the mask from charging. Using an adapter with voltage above 5 V may damage the product.

— This device features overvoltage protection (for input voltages between 6V and 20V).

— Never use the product while charging.

Disposal

1. When disposing of the device, remove the battery first.
2. Ensure the device is powered off before battery removal.
3. After extraction, dispose of the battery according to local regulations.

Warranty

The Rejuvenating Phototherapy Mask LED MASK 2.0 by WAU comes with a 1-year warranty. During the warranty period, the product will be repaired or replaced if any manufacturing defects are found. This warranty does not cover faults caused by improper use or transportation, or products that have been self-repaired.

Frequently Asked Questions

Q: Is this mask waterproof?

A: No. Only wipe gently with water or ethanol. Do not wash or submerge.

Q: How long can the mask operate on a full charge?

A: In mixed (white) light mode, one charge provides 5-6 sessions of 15 minutes each. In red/blue/green light modes, the mask can be used for 8-9 sessions of 15 minutes.

Q: My skin became dry after using the mask. How can I fix this?

A: Since the light wavelengths penetrate to the dermal layer and stimulate cell metabolism, enhancing the skin's absorption capacity, your skin requires more hydration and nourishment. Therefore, after using the LED MASK 2.0, we recommend replenishing moisture: provide adequate water and nutrients to your skin.

Q: Does the mask heat up during use?

A: When using red, blue, or green light modes, the mask temperature remains stable. In mixed (white) light mode, the mask may warm up to $38\pm 1\text{ }^{\circ}\text{C}$ ($100,4\pm 1,8\text{ }^{\circ}\text{F}$).

Q: Are the results visible immediately after the first use?

A: The effects after one or two applications will be less noticeable because skin metabolism follows a 28-day cycle. Therefore, visible results typically appear after 1 month of consistent use.

Q: What is the recommended session duration for the mask?

A: Limit phototherapy sessions to 15 minutes maximum. Exceeding this duration may cause skin damage. Please monitor your usage time carefully.

Q: Can I bend or fold the mask?

A: No, as this may damage its internal circuitry.

LED MASK 2.0

REJUVENATING PHOTOTHERAPY MASK

Model No.: **K9C-1KB**

Material: **ABS-Kunststoff, Silikon**

Dimensions (H×W×D): **190x679x5 mm**

Weight: **380 g**

Voltage: **3,7 V**

Power: **3,7 W**

Manufactured in: **PRC**

Authorized manufacturer representative: WAU Global Limited
63A College road, Northcote, Auckland, 0627, New Zealand
info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com