wau

LED MASK

PHOTOTHERAPY FACE MASK

USER GUIDE

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Warranty 31

Technical Specifications


Parameter	Value
Model	FM-22
Materials	ABS plastic, silicone
Dimensions (H×W×D):	198x383x33,5 mm
Weight	216,5 g
Power Source	Lithium-ion polymer battery
Battery Capacity	1500 mAh
Operating Voltage	3,7 V
Power Consumption	3,7 W
Charging Time	3-4 house
Operating Time	40-45 minutes
Modes	Red / Blue / Yellow / Mixed (Purple) light
Number of LEDs	100 lamps
Irridation Intensity	• Red light — 2,4041 mW/cm² • Blue light — 8,6094 mW/cm² • Yellow light — 1,2704 mW/cm² • Mixed (Purple) light — 7,4833 mW/cm²


Purpose

The LED MASK is a home-use phototherapy mask from the brand WAU — a unique device that delivers positive effects to facial skin through the use of LEDs emitting different wavelengths: red, blue, yellow, and mixed light.

Depending on the wavelength, different colors of light penetrate various layers of the skin, providing beneficial effects that make the skin appear more even, hydrated, and smooth. Additionally, LED therapy can help accelerate metabolism, slow down skin aging, reduce sebum production, effectively disinfect, combat acne and inflammation, reduce swelling, and improve complexion.

Device Components

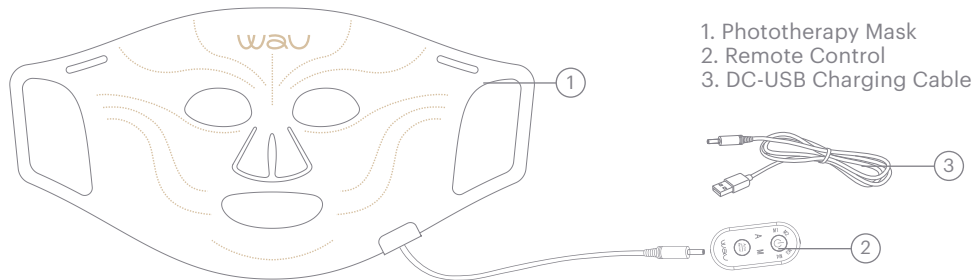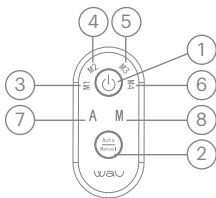Remote Control Labels


1. Power On/Off/Mode Button
2. Manual/Automatic Mode Button
3. Red Light
4. Blue Light
5. Yellow Light
6. Mixed (Purple) Light
7. Automatic Mode Indicator
8. Manual Mode Indicator

Package Contents

Phototherapy Mask ×1, Remote Control ×1, DC-USB Charging Cable ×1, Mounting Straps ×4, Eye Protectors ×2, User Manual ×1.

Operating Procedure

1. Assemble the mask by securing the straps according to the illustration.

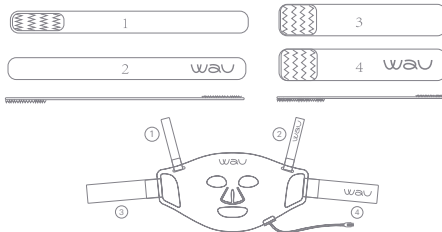2. On the remote control, press and hold the Power On/Off button for 1 second to turn on the mask. Upon activation, the red light function and manual mode will be enabled. An indicator with the letter M (Manual mode) will light up next to the mode switch button. If the mask requires charging, please refer to the Charging section of the manual. Note that the mask does not operate while charging.
3. Press the Power button again to switch the mask to blue light mode.
4. Press the Power button a third time to activate yellow light mode.
5. Press the Power button a fourth time to enable the mixed phototherapy mode.
6. In manual mode, the mask operates for 10 minutes.
7. To switch the mask to automatic mode (automatic light switching), press the mode switch button. In automatic mode, the mask cycles through: Red light for 5 minutes, Blue light for 3 minutes, Yellow light for 4 minutes, Mixed (Purple) light for 3 minutes
8. To turn off the mask, press the Power On/Off button.

Usage Method

[Important Note] Before first use, it is strongly recommended to perform a sensitivity test:

- Place the device on the inner side of your arm;
- Turn on the mask and run a full 10-minute phototherapy session;
- If you experience any discomfort or burning sensation, discontinue use immediately;
- After the session, wait at least 6 hours. If redness and/or itching occurs on the tested area, the LED MASK is not suitable for your use.

1. Wash and thoroughly dry your hands and face.
2. Put on the mask, ensuring that the eye protectors fit tightly around your eyes.
3. Select either manual or automatic mode as desired.
4. In manual mode, choose your preferred light therapy color.

Note: When switching light modes in manual operation, the 10-minute smart timer resets. However, it is strongly advised not to exceed 15 minutes of use per session.

5. Enjoy your phototherapy session in the comfort of your home.
6. After the session, it is recommended to apply a moisturizing serum or cream, as the skin's absorption ability is enhanced post-treatment.
7. For visible results, use the mask regularly — 3 to 4 times per week for at least 4 weeks.
8. Once the desired effect is achieved, maintenance sessions can be performed 1 to 2 times per week.

Mode Features

Red Light (630 nm)

Having the longest wavelength, red light penetrates the skin deeper than the other presented colors. It acts on fibroblasts, stimulating collagen production and thereby promoting skin regeneration. As a result of red light therapy, wound healing is accelerated, the skin becomes firmer, and visible signs of aging such as wrinkles are reduced.

Blue Light (490 nm)

Blue light has the shortest wavelength among all visible spectrum colors. It works on the level of the epidermis, providing antimicrobial effects. Blue light is used to treat acne by destroying the bacteria that cause it.

Yellow Light (590 nm)

A type of light that penetrates the skin deeply enough to stimulate collagen production, helping to smooth wrinkles. Yellow light improves lymphatic flow and drainage, which helps reduce puffiness. In addition, it enhances skin tone, lightens spots, and speeds up the reduction of pigmentation.

Mixed (Purple) Light

The combination of red and blue light merges their functionality. Blue light addresses skin conditions related to the outer layer of the skin — the epidermis. On the other hand, red light significantly influences sebum production. Thus, mixed light is the most effective for treating inflammation.

Precautions

- The LED MASK is a device intended for home use. It may only be used in accordance with the method described in this user manual. It is prohibited to use it for other purposes, as this may result in product damage or harm to the body.
- Check the device for any malfunctions before each use.
- If the mask is damaged, stop using it and contact the authorized consumer service center to avoid any danger. In cases of self-repair, disassembly, or improper use, the authorized consumer service center is released from warranty obligations.
- Do not use the mask while it is charging.
- Do not touch the mask, charging cable, or USB adapter with wet hands during charging.
- Do not allow water to get on the surface or inside the device.
- Stop using immediately if you experience discomfort or any side effects (skin redness, itching, allergic reactions) as a result of using the mask.
- The device is not intended for use by infants, people with epilepsy, individuals with thyroid diseases, or those taking thyroid medication.
- The mask must not be used on damaged, injured skin or post-operative wounds.
- Do not use the device near open heat sources.

Storage

1. Clean the surface of the LED MASK after each use.
2. Store the device in a dry place at room temperature.
3. Avoid storing the device in dusty areas, under direct sunlight, or in a humid environment.
4. This mask contains no user-serviceable parts.

Cleaning & Maintenance

1. Disconnect the device from the power source during the cleaning process.
2. Do not spray water directly onto the device when cleaning.
3. Wipe the surface of the mask twice (three minutes each time) using a 70–80% ethanol solution. Then wipe with a damp or dry cloth.
4. Do not use flammable substances such as benzene, solvent, butyl acetate, etc., for cleaning this product.
5. Do not use other cleaning agents, as this may damage the mask.
6. Do not rinse the device directly under running water during cleaning.

Charging

- Insert the DC charging cable connector into the corresponding port on the remote control, and plug the other end with the USB connector into the power adapter.
- Ensure that the power adapter's output voltage matches the voltage specified on the device.
- At the start of charging, the M1 indicator on the remote control will light up. It will blink until the charge reaches 25%, then the M1 indicator will stay lit continuously, and the M2 indicator will start blinking. When the charge reaches 50%, indicators M1 and M2 will remain lit steadily, and indicator M3 will blink. At 75% charge, indicators M1, M2, and M3 will stay lit continuously, and indicator M4 will blink. After charging is complete, all indicators will remain steadily lit.
- Using a power adapter with an output voltage below 5 V will not charge the device. Using a power adapter with an output voltage above 5 V may damage the remote control.
- If the remote control battery voltage drops below 3,2 V, the remote control indicator will start blinking, and after 5 flashes, the device will automatically shut down.
- The mask is protected against overvoltage (protection activates between 6 V and 20 V). Use of the device during charging is prohibited.

Disposal

1. When disposing of the device, the battery must be removed.
2. The device must be turned off before removing the battery.
3. After removing the battery, dispose of it in accordance with local regulations.

Warranty

The LED MASK phototherapy mask by WAU is covered by a 1-year warranty. During the warranty period, the product is subject to replacement or repair in case of technical defects. This warranty does not cover malfunctions caused by improper use or transportation, or products that have been repaired independently.

LED MASK

PHOTOTHERAPY MASK

Model: **FM-22**

Material: **ABS plastic, silicone**

Dimensions (H×W×D): **198 × 383 × 33,5 mm**

Weight: **216,5 g**

Voltage: **3,7 V**

Power: **3,7 W**

Made in **PRC**

Authorized manufacturer representative: WAU Global Limited 63A College road,
Northcote, Auckland, 0627, New Zealand
info@wauglobal.com

Important Notice: to ensure continuous product improvement, the manufacturer reserves the right to modify design and technical specifications. The latest version of this manual is available at: wauglobal.com