wau

VIA Gold
LIFTING MASSAGER

USER GUIDE

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RoHS

Table of Contents

Technical Specifications.....	18
Purpose	19
Device Components.....	19
Operation Procedure.....	20
Usage Method	21
Cleaning & Maintenance.....	25
Storage.....	25
Charging.....	25
Precautions.....	26
Warranty.....	29


Technical Specifications

Parameter	Value
Model	MR-1821
Materials	ABS Plastic, PC Plastic
Dimensions (HxWxD)	173x55x45 mm
Weight	208 g
Power Source	Lithium-ion Battery
Battery Capacity	1000 mAh
Operating Voltage	3,7 V
Power	3 W
Charging Time	4 hours
Operating Time (full charge)	SONIC Mode 150 minutes RF Mode 40 minutes EMS Mode 2700 minutes
Smart-Timer	SONIC Mode 3 minutes RF Mode 5 minutes EMS Mode 6 minutes
Water Resistance Rating	IPX6

Purpose

VIA Gold is an effective skin rejuvenation massager that combines the best hardware technologies: RF, EMS, and Sonic. It is specially designed for performing comprehensive facial and neck care procedures. The SONIC mode completes the cleansing stage, RF promotes collagen regeneration, restoring skin firmness and elasticity, while EMS stimulates metabolic processes to lift the facial contour and improve muscle tone.

Device Components

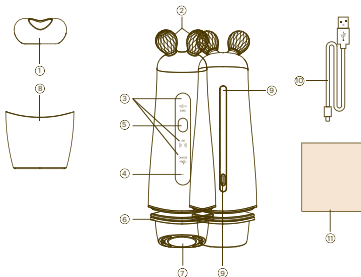1. Transparent Protective Cap
2. V-Shaped Face Rollers
3. Mode Indicators
4. Battery Indicator
5. Power/ Mode Button
6. Magnetic Ring
7. Applicator
8. Magnetic Base
9. Induction Strip
10. USB-Type-C Charging Cable
11. Cleaning Cloth

Package Contents

Device x1, Magnetic Base x1, Protective Cap x1, Magnetic Ring x1, USB-Type-C Charging Cable x1, Cleaning Cloth x1, Storage Pouch x1, User Guide x1.

Operating Procedure

1. Remove the protective cover and cap from the device.
2. Turn on the massager by pressing and holding the Power/ Mode button.
3. The device offers three modes: 1 — SONIC; 2 — RF; 3 — EMS. Switch between the modes by briefly pressing the Power/ Mode button. A short beep will accompany each mode change.

4. In EMS mode, you can adjust the intensity level by pressing the Power/ Mode button again. The brightness of the indicator light reflects the intensity: the brighter the light, the stronger the current.
5. The device features a smart timer: 3 minutes for SONIC mode, 5 minutes for RF mode, and 6 minutes for EMS mode. Once the set time is reached, the device will turn off automatically.
6. To turn off the device manually, press and hold the Power/ Mode button. A long beep will sound both during automatic and manual shutdown.
7. The VIA Gold massager is equipped with a memory function: the device will automatically start in the last-used mode upon the next activation.

Usage Method

SONIC Mode

1. Remove makeup and cleanse your face using a facial cleanser.

2. Apply toner to a cotton pad or use a cleansing pad.
3. Secure the pad onto the applicator using the magnetic ring.

4. Turn on the device and select the SONIC mode.

5. Gently apply the applicator to your face and begin the cleansing procedure, moving along the massage lines from the bottom of the chin upward toward the hairline on the forehead. For best results, cleanse each area multiple times, paying special attention to the T-zone.

6. After 3 minutes of operating, the device will automatically turn off.

7. Remove the pad and clean the applicator with the cleaning cloth.

8. Rinse your face with water.

RF Mode

1. Apply a water-based skincare product to your face. Ensure the formula does not contain aggressive acids or oils.

2. Turn on the device and switch to RF mode.

3. Slowly move the applicator from the center of the face outward, maintaining a consistent speed. Do not hold the applicator in one spot for too long.

4. After 5 minutes, the device will automatically switch off.

5. Use the cleaning cloth to remove any remaining product from the applicator.

EMS Mode

1. Apply a water-based, gel-textured skincare product to your face. This will improve electrical conductivity for optimal results.

2. Turn on the device and switch to EMS mode.

3. Select the EMS intensity level — low or high — by pressing the Power/ Mode button again.

4. Place the rollers along the facial contours and gently massage upwards, slightly lifting the skin. Perform the procedure 3–5 times, moving along the massage lines from the center of the face outward. Avoid the forehead and eye areas.

5. After 6 minutes, the device will automatically turn off.

6. Use the cleaning cloth to remove any remaining product from the rollers.

Cleaning & Maintenance

After each use, rinse the applicator and V-shaped rollers with clean water, then dry thoroughly with a clean cloth. Do not use alcohol or corrosive liquids to clean the magnetic ring.

Storage

- Store the device in a dry place, protected from dust.
- Keep the device away from high temperatures.
- Keep the device out of reach of children.
- Protect the device from exposure to corrosive liquids.
- Do not immerse the device body in water.

Charging

- When the battery charge falls below 10%, the white battery indicator will flash 3 times every 3 seconds. Before the battery is completely drained, a double beep sounds, the indicator starts flashing rapidly, and the device turns off.
- At 0% battery charge, the device emits a triple beep when turned on. The battery indicator flashes for 5 seconds, then the device powers off. It cannot be turned on again until charged.

- Connect the Type-C charging cable to the device's charging port, and plug the USB end into an adapter.
- During charging, the white battery indicator will slowly turn on and off.
- When fully charged (100%), the white battery indicator stays steadily lit.

Precautions

Do not use the device:

- If you have undergone surgery within the past month, or if you have severe sunburn on your skin;
- If your skin is inflamed, scratched, damaged, or if you have an allergy to electromagnetic waves;
- By children or individuals with blood disorders, tuberculosis, internal organ diseases, malignant tumors, high blood pressure, high body temperature, or facial nerve neuralgia;
- In EMS mode around the eyes, mouth, thyroid area, or on areas with scars from plastic surgery or injuries;
- On the same area for an extended period.

Before use:

- Remove all metal items such as watches, glasses, etc.

During Use

— Do not wear clothing made of synthetic fibers while using the device.

With Caution

— If you have facial silicone implants or dental metal crowns, pins, or braces, use the device carefully. Pay attention to your own sensations and skin sensitivity. If the device causes significant discomfort, discontinue use immediately.

Important

- Do not leave the device unattended while it is operating.
- Do not use the device in high-temperature environments.
- Do not disassemble the device yourself.
- Do not turn on the device while it is charging.
- After dropping or striking the device, check for damage before continuing to use it.
- Regularly inspect the plug and power cable to avoid the risk of electric shock.

Malfunctions and ways to eliminate them

Malfunction	Reasons
The device does not work after switching on	1. The power indicator does not light (battery or device circuit failure) 2. The power indicator flashes quickly (battery is low)
How to fix	
1. Charge the device for 4 hours, then turn it on again. If the power indicator does not light, contact the after-sales service and return the device to the organization authorized by the manufacturer to accept consumer claims. 2. Charge for 4 hours, then turn it off again.	

Malfunction	Reasons
The device charges poorly	1. Damaged charging cable 2. Damaged USB port 3. Battery malfunction
How to fix	
1. Make sure the charging cable is properly connected, the charging indicator is blinking, and charging is in progress. 2. If nothing changes after charging, contact an organization authorized by the manufacturer to handle consumer claims.	

Malfunction	Reasons
The device shuts down unexpectedly while in use	1. The battery is discharged 2. Other issues have occurred
How to fix	
1. Charge for 4 hours, then turn the device off. 2. Contact an organization authorized by the manufacturer to handle consumer claims.	

Recommendation: If the device has not been used for an extended period (6 months or more), it should be fully charged before use.

Warranty

The VIA Gold Lifting Massager by WAU is covered by a 1-year warranty. During the warranty period, if technical defects are found, the product is subject to replacement or repair.

This warranty does not cover malfunctions caused by improper use, transportation damage, or products that have been repaired by unauthorized persons.

VIA Gold

LIFTING MASSAGER

Model: **MR-1821**

Material: **ABS plastic, PC plastic**

Dimensions (HxWxD): **173 x 55 x 45 mm**

Weight: **208 g**

Voltage: **3,7 V**

Power: **3 W**

Made in: **PRC**

Authorized manufacturer representative: I. E. Umantsev A. N.
Leningradsky Prospekt 37, Office 14-14, Moscow, 125167, Russia.
info@wauglobal.com

Important Notice: to ensure continuous product improvement,
the manufacturer reserves the right to modify design and
technical specifications. The latest version of this manual is
available at: wauglobal.com