wau

ELAX 2.0
EYE RELAX MASK

USER GUIDE

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Table of Contents

Technical Specifications	22
Purpose	23
Device Components	23
Operating Procedure	24
Usage Method	25
Cleaning & Maintenance.....	26
Storage	26
Charging	26
Precautions	27
Warranty.....	28


Technical Specifications

Parameter	Value
Model	LX-2
Materials	Modal, Lycra
Dimensions (HxWxD)	108x260 mm
Weight	111 g
Power Source	Lithium battery
Battery Capacity	1000 mAh
Operating Voltage	5 V
Power	10 W
Charging Time	1 hour
Operating Time (from a full charge)	75 minutes (high temperature mode) 60 minutes (medium temperature mode) 45 minutes (low temperature mode)
Smart-Timer	15 minutes
Heating Temperature Levels	Low: 40-45 °C Medium: 45-50 °C High: 50-55 °C

Purpose

The ELAX massage and relaxation mask by WAU offers comfortable use without pressure or heaviness. Three adjustable heating modes and a moderate vibration level help relax muscles, improve blood circulation, and relieve eye strain. The mask comes with a cold compress to reduce puffiness. The 3D design ensures no pressure on the eyes or eyelashes, provides a snug fit to the face, and blocks out light, allowing the device to be used as a sleep mask as well.

Device Components

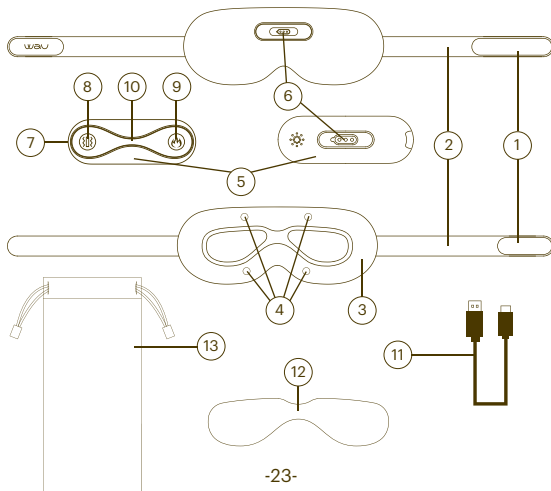
1. Velcro fasteners
2. Elastic strap
3. 3D mask base with memory foam function
4. Vibration motors
5. Control panel
6. Control panel holder
7. Charging port
8. Vibration mode on/off button
9. Heating mode on/off and adjustment button
10. Battery level indicator
11. USB-C charging cable
12. Cold compress
13. Storage pouch

Package Contents

Mask x1, cold compress x1, USB Type-C charging cable x1, storage pouch x1, user manual x1.

Operating Procedure

1. Connect the control panel to the mask.
Note: The battery indicator will light up white if the connection between the mask and control panel is weak.
2. Press and hold (for 2 seconds) either the vibration power button or the heating mode button to turn on the device. The battery level indicator will light up.
3. The vibration and heating modes operate independently. Select the desired mode.

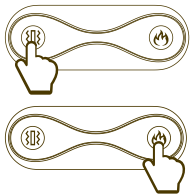4. The vibration mode has three intensity levels. The massager starts operating at the low vibration level, emitting one sound signal. When switching levels, two sound signals indicate the medium vibration intensity, and three sound signals indicate the maximum vibration intensity.

5. The heating mode has three temperature levels. The device starts at the low heat level, accompanied by one beep. Switching to the medium heat level triggers two beeps, and the high heat level triggers three beeps.

6. The device features a smart timer that automatically shuts off after 15 minutes of operation. The vibration and heating modes have independent timers that start separately for each mode. Changing the temperature level does not reset the smart timer.

7. To turn off the device, press and hold the vibration or heating power button for 2 seconds.

Usage Method

1. Connect the control panel to the mask.
2. Put the mask over your eyes and turn it on.
3. Activate the vibration mode, heating mode, or both simultaneously.
4. Select the desired temperature level for the heating mode.
5. Enjoy a 15-minute massage session.
6. If you wish to use the cold compress, place the insert in the refrigerator for 10 minutes beforehand. Then insert it into the mask without turning the device on.

Cleaning & Maintenance

- If you need to clean the mask, remove the control panel and wipe the mask with a damp cloth or one moistened with water.
- Do not scrub the mask vigorously.
- After cleaning, allow the mask to air dry naturally. Do not use dryers, hairdryers, or other heating devices.

Do not wash

Do not iron

Do not bleach

Dry cleaning
prohibited

Do not wring
or tumble dry

Storage

- Store the device in a dry place, protected from dust.
- Keep the device away from high temperatures.
- Avoid exposing the device to direct sunlight.
- Store the device out of reach of children.
- If the device is not used for an extended period, turn it on and charge it once every 2 months.

Charging

- When the battery level is above 66%, the battery indicator lights up green.
- When the battery level is between 33% and 66%, the indicator lights up orange.

- When the battery level is below 33%, the indicator turns red. The device will continue to operate, and the red indicator will blink during use.
- Connect the Type-C charging cable to the device's charging port. Plug the USB end of the cable into a power adapter.

- During charging, the battery indicator blinks red.
- The device does not operate while charging.
- Once fully charged (100%), the battery indicator lights up green.

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- Keep the product out of reach of children and infants.
- This product is intended for users aged 12 and above.
- Keep the control panel, charging cable, and charging port dry.
- Do not spray chemical liquids such as perfumes or deodorants on the surface of the mask.
- Do not bend, pull, or tie the eye mask or charging cable to avoid damage.
- Individuals with eye conditions, allergies, or sensitivity to heat compresses are advised to consult a doctor before use.

Warranty

The ELAX Eye Relaxation Mask by WAU is covered by a 1-year warranty. During the warranty period, the device is eligible for replacement or repair in the event of any technical defects. This warranty does not apply to malfunctions caused by improper use, transportation damage, or unauthorized repairs performed by the user.

ELAX 2.0

EYE RELAX MASK

Model: **LX-2**
Materials: **Modal, Lycra**
Dimensions (H×W×D): **108×260 mm**
Weight: **111 g**
Voltage: **5 V**
Power: **10 W**
Made in: **PRC**

Authorized manufacturer representative: PERFECT SQUARE
LIMITED 54/2 Fairclough Rd, Beach Haven, Auckland,
Auckland, 0626, New Zealand
e-mail: info@wauglobal.com

Important Notice: to ensure continuous product improvement,
the manufacturer reserves the right to modify design and
technical specifications. The latest version of this manual is
available at: wauglobal.com

RoHS